

Breast Enhancer Operating Instructions

Advice of experts:

Do be cautious to select right breast enhancing products

Today there are a lot of products named breast cream, breast paste in the shelf of breast enhancing series in the market, a fraction of them contain estrogen especially diethylstilbestrol, applying it maybe enlarge the breasts to some extent, but it will remain a short time only, further, once you stop using the medicine, the breasts immediately return to original form, in addition, such medicines may result in pigment sediment, black freckle, abnormal menstruation and other side effects, For any teen girls of normal menstruation, they have sufficient estrogen in their body, but if we considerably add estrogen so as to enhance the breast, it will restrain the secretion of internal estrogen, and thus restrain the growth of breast, it is just the opposite to what one wishes.

Taking orally the medicine for breast enhancing may be worsened, because it aims to gain larger breast by stimulating the secretion of estrogen in body, but while the breast is extending, the waist is also extending. Additionally, taking a large dosage of breast enhancing medicine damages your liver, kidney and stomach, long-term taking it may seriously injure your health.

Three features of this product

1. As a micro-computer high-tech product, it is one of the lightest and easiest health apparatus in the world.
2. Its controller employs high-tech chip, it is safe, reliable and endurable without irritation or rebound.
3. It is a breast care expert without any side effect.

An amazing move generates a new age of breast enhancing. Expert tips

Psychologically the women with small and flat breasts mostly feel self-abased, they often suffer other mental disease or encounter early aging, early turn of life and other symptoms.

Sexologically, the size of breast is closely related to the level of estrogen, most of the women of flat breast are undersexed, really, for married women, this will seriously affect the quality of life.

Medically, during lactation period, the maldeveloped breasts may encounter disability in lactic secretion, it will not only disturb lactation, but also cause an array of diseases including parotitis.

Application:

1. Flat breast left in maldeveloped adolescence or caused by innate disability, or innutriton.
2. Shrink, flat, slack, droppy, collapse breast after swelling during pregnancy or lactation.
3. Seriously shrink or droppy breast because of poor secretion of prolactin or progesterone with age increasing, or under effect of gravity, thus cause the skin to lost bounce, breast muscle to become slack.
4. The woman at age of 18-45 and the growth of breast fails to meet their expectation.

Working mechanism

Dredging- By introducing advanced high-tech micro-biological electricity and using the current in different modes, dredge the channels of breast, enhance micro-circulation of breast, improve vitality of cells, It clearly prevents any induration or other pathology in breast.



Actuating- The unique micro-electricity and its different strength modes stimulate the pituitary of hypothalamus to release gonadotropic hormone, also stimulates ovary to secrete estrogen progesterin, enhance the responsiveness of breasts to internal hormone, resulted in better growth of mammary gland and breast lobules, and make the breast vein quickly extend.



and ramify, so as to enable breasts to expand and swell rapidly.


Slimming- You can repetitively press ON/▲ button to increase strength, the moment you will feel a little irritated in the massaging parts, then you can readjust the strength to desired level. Using the device 5-15 minutes every day can efficiently enhance tenacity of breast suspensorium and its net structure, so as to raise the droopy breast, make your breast tough and in good shape.



Maintaining- Long-term using this product can successfully adjust the hormone secretion, maintain the sensible responsiveness of your body to hormone, keep steady the nutrition of tissues and cells, make your breast plump and elastic lastingly.



Diagram of host

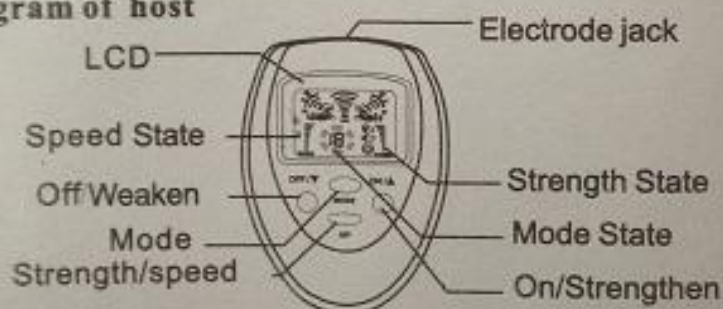

Preparations

1. Load-in batteries

- 1) Slide down the battery cartridge to open the cover.
- 2) Place batteries to the cartridge with label side toward you.
- 3) Close the battery cover.


2. Connections

- 1) Set the connecting cables and electrode pads as (A) and (B), You can use one or two pads at the same time with all cables, Pins to receptacles marked Breast Enhancer.


2) Place electrode pads on fatty part of body as (D)

For Breast Developing


1. Adjust the silicon cups to an equal size. Apply breast enhancing oil to the black parts of silicon cup, and place the silicon cup between breast and bra as (D).
2. Insert the connecting pins into silicon cup and main body as (B), (C).
3. Fix the set as (E).
4. To strengthen stimulation performances, different electricity current and pause duration combination can be adjusted as you like to most comfortable situation.

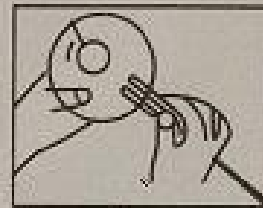

A

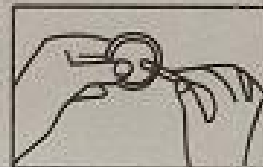

B

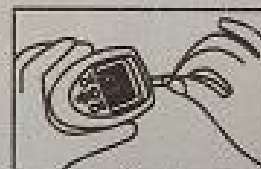

C

BOSOM DEVELOPER:

1. Recreation for every woman with sexier shape.
2. Broadening under-sized ones.



E

3. Toughening droopy breasts.
4. Promoting complete development of breasts.
5. Mothers just after the accouchement to increase lactation.

NOTE & PRECAUTIONS:

1. When finished using always turn power knob off first.
2. Clean skin thoroughly before use, Clean any dirt, lotion, oil and make-up left in skin. Apply only this product to unbroken skin.
3. After finished using return electrode pads and silicon cups in the box at normal temperature.
4. Do not use under following physical conditions:
 - Skin disease
 - Abnormal blood pressure
 - Pregnant
 - Brain hemorrhage
 - Heart disease
 - Malignant tumor
5. Immediately stop using when you feel any physical discomfort.