



MS-3351B

Instructions

Warning

Dear users:

Thanks for choosing our company's product, in order to fully make use of the products, we sincerely suggest follow:

- ◆◆Please read the instruction in detail and keep it well for further reading and consulting.
- ◆◆Please follow the instruction to fit and operate the machine correctly.
- ◆◆Please do not remove and change any accessories of the machine.
- ◆◆It forbids others to open and remove the mainframe except the authorized people in the company.
- ◆◆It forbids others to change and maintain the machine except the professional people.
- ◆◆Please contact us if there is any trouble during operation.

After Sale service

◆We offer ONE YEAR warranty. Within warranty period, all exchanges or replacement due to defects are free. (Buyer is responsible for shipping charges).

If the machine is over warranty, please also do not hesitate to contact us. we will assure you of our best service all the time.

◆How to contact us?

Phone:(+86) 189 -487-55929

Tel: 0755-86368500

Email: sales@mychway.com

This product does not belong to medical, please don't use for disease treatment!

Chapter 1: Safety Warnings

I. Electricity & Machinery Safety

1. The machine is either 110V And 220V. So whatever the voltage u can use before starting the machine.
2. Please keep the voltage at stable condition.

II. Contraindication Safety

1. People with heart disease, high blood pressure, or configured cardiac Pacemaker.
2. People with acute inflammation, asthma, deep vein thrombosis, thyroncus, cancer.
3. People with hemorrhagic diseases, trauma or who is bleeding.
4. Pregnant women
5. Medical plastic artificial or metal parts inside the body such as metal teeth or silicone breast implants.
6. People with abnormal immune system.
7. People with numb or insensitive to heat.
8. If you are feeling unwell, Prohibited in the ears, nose, eye throat and do not use this device in advance and tell your doctor.
9. Children under 12 years old are not allowed.

III. Packing Including:

- 1 x Main Machine
- 1 x Unoisetion Caviation 40K Probe
- 1 x Quadropole 3D RF Probe for face/body
- 1 x Bipolar 3D RF Probe for face and eyes
- 1 x Holder
- 1 x Power Supply Cord

IV. Treatment Safety

1. Treatment Time is less than 30mins/Body Part/Each Treatment.
2. The energy setting is from low to high gradually which should be subject to skin

comfortableness for different people.

3. The patient had better wear the protective eye glasses.

4. If you cannot tolerant heat or uncomfortable, please stop the treatment in advance and turn to your doctor for help.

Chapter 2: Installment and Use

I .Working Condition

Using at normal temperatures and pressure. Please do not stay in too cold, too hot, too dry and too wet.

II. Technical Parameters

Power supply Input:100V-240V

Power :110W

Unoisetion Probe:

Power supply output:150V

Frequency:40KHz

Power:25W

Quadrupole 3D SMART RF head for Arm/face/body

Power supply output:66V

Frequeuncy:3MHz

Power:30W

Bipolar 3D SMART RF for face and eyes

Power supply output:66V

Frequeuncy:3MHz

Power:30W

III. Installation



Bipolar RF handle



Quadrupole RF handle



40K Unoisetion cavitation handle




Power supply cord

Chapter 3: Operation


1. Connect the power supply.
2. Turn on the red button of power switch.



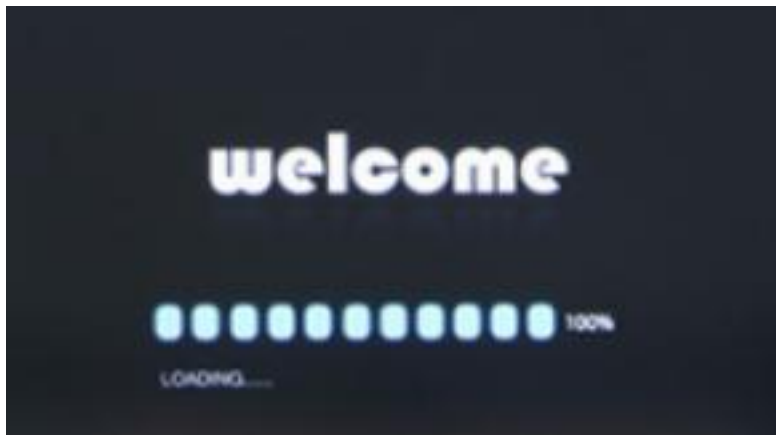

3. When turn on machine, you will see screen show like this:




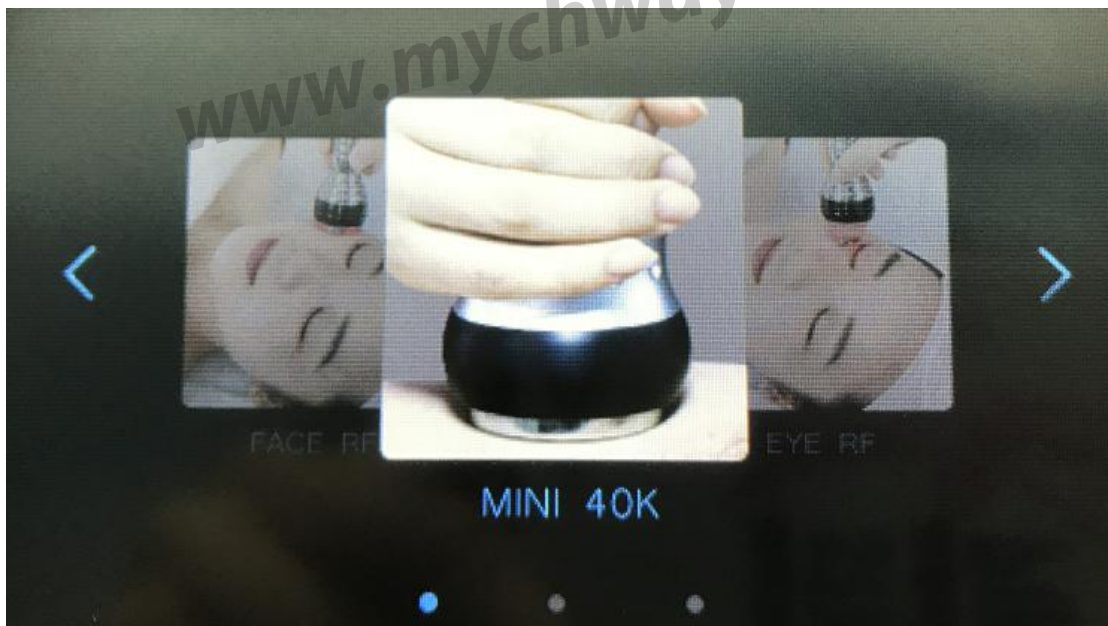










MINI 40K

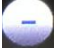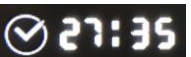


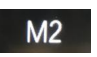

1. Click 40k picture to enter cavitation function, or click  and  to choose face rf or eye rf function.





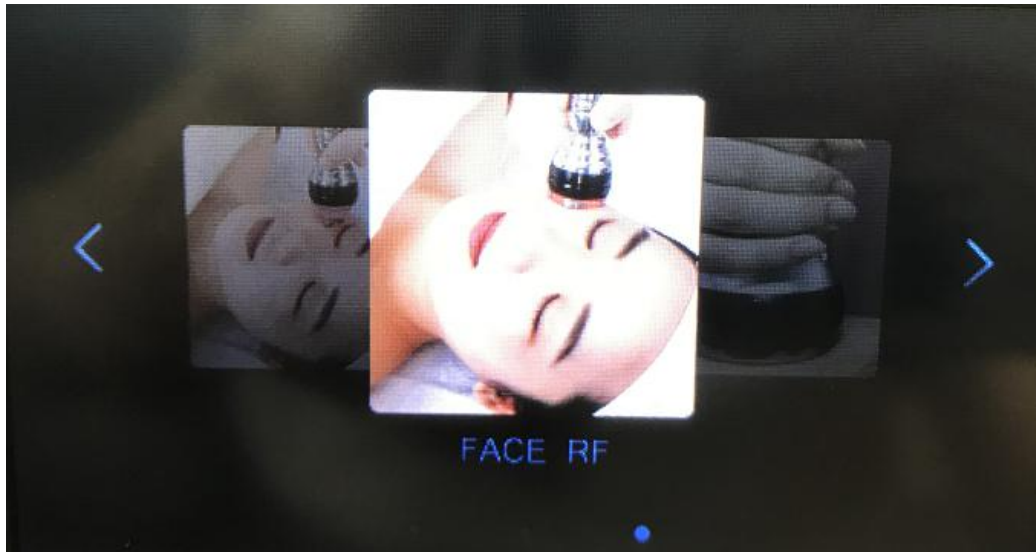








2.  Click here, and press  or  to plus/ minus the energy needed.
3.  click time setting, and press  or  to plus/ minus the time treatment needed. Default time is 30mins.
4.  Press here, and press  or  to choose the mode you need, there are M1 M2, 2 modes. M1 means continuous working mode , M2 means pause working mode, please chose the most comfortable mode upon your feeling. Default mode is M1.
5.  turn on /off the function.
6.  Back to home page.
7.  Reset your setting.

Face RF




1. Click face rf picture to enter into the function, or click  and  to choose 40k or eye rf function.





2.  Click here, and press  or  to plus/ minus the energy needed.
3.  **30:00** click time setting, and press  or  to plus/ minus the time treatment needed. Default time is 30mins.
4.  click here, and press  or  to turn on / off the led light on face rf.
5. **NOR PRO** click here, and press  or  to choose the mode of face rf, there are NOR or PRO, two modes.

Nor: means normal mode for this RF working model, temperature will reach the one setted at once.


Pro: professional mode, temperature will reach to the one u setted slowly.



6.  turn on /off the function.

8.  Back to home page.

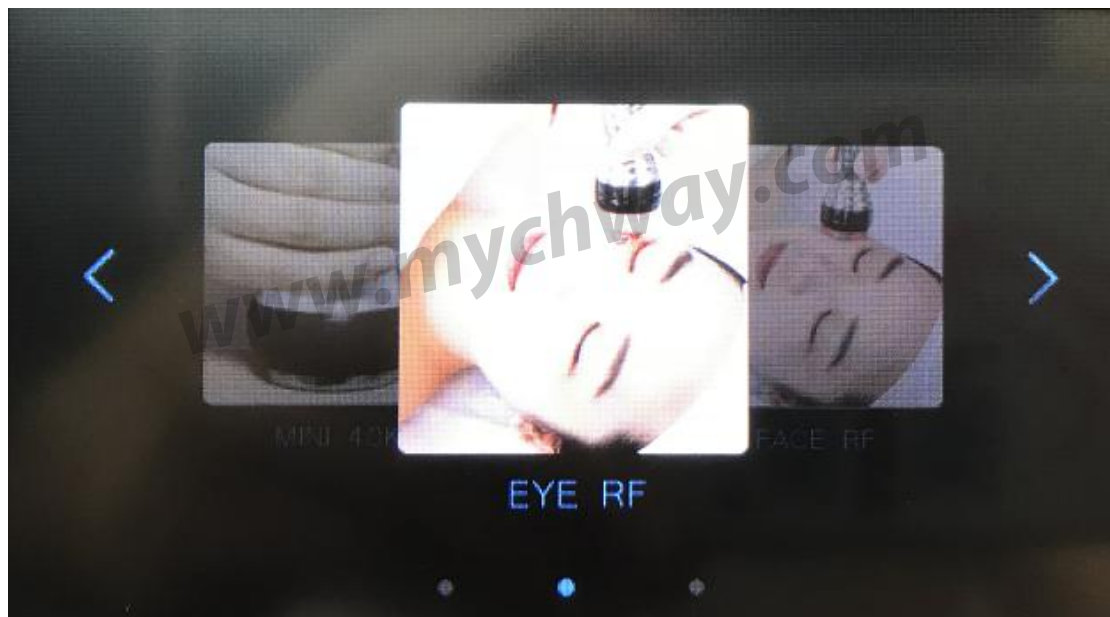


9.  Reset your setting.




EYE RF for Eye






1. Click eye rf picture to enter into the function, or click  and  to choose 40k or face rf function.



7.  Click here, and press  or  to plus/ minus the energy needed.


8.  click time setting, and press  or  to plus/ minus the time treatment needed. Default time is 30mins.

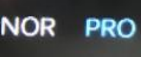
9.  click here, and press  or  to turn on / off the led light on eye rf.


10.  click here, and press  or  to choose the mode of face rf, there are NOR or PRO, two modes.

Nor: means normal mode for this RF working model, temperature will reach the one setted at once.

Pro: professional mode, temperature will reach to the one u setted slowly.

11.  turn on /off the function.

10.  Back to home page.

11.  Reset your setting.

Chapter 4: Treatment time and operating

- 1) Abdomen, thigh, buttock, arms slimming

Please clean the skin and then put the ultrasound gel on the treatment area, then turn on the machine, set the time and energy. Please set the energy from low to high to get the best comfortable temperature.

1. Please use the unoisetion 40k cavitation head 25-30min

Work time: 2-3 times/ week

1. Face, eyes around, neck skin lifting:

Please clean the skin and then put the rf gel on the treatment area, then turn on the machine , set the time and energy, pls set the energy from low to high to get the best comfortable temperature, then follow this steps:

A: Use face radio frequency head for face for 25-30 min


B: Use eye radio frequency head for eye around for 15-20min

Work time: 2-3 times/ week

OPERATION

FACE

Treated time: 30 minutes

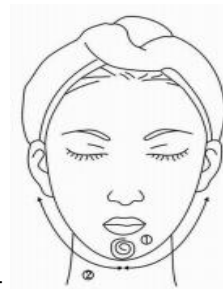

1. Massage the mandible center in circles to produce deep heat.

2. Massage lower jaw in lines.



3. Massage the triangle zone of both sides of the face in circles.

4. Massage from jaw to angulus oris to ear in lines.

THIGH

Treated time: 30 minutes



1. From down to upper, pushing to the groin to dredge the lymph.

2. From down to upper, circling push by anticlockwise can help decomposing



fatness.



3. Also can push by come-and-go to decompose fatness.

4. Pull from the knee and muscle texture to upper, can improve the curve.

Treated time: 20 to 30 minutes

1. Along the lymph direction, pulling to the waist.



2. From upper to down, by anticlockwise gesture, pull come-and-go.



3. Along the muscle of arm, pull up to lift and tight the muscle.


BACK

Treated time: 20 to 30 minutes

www.mychway.com



1. First by come-and-go, pushing the bladder nerve 2 to 3 times.

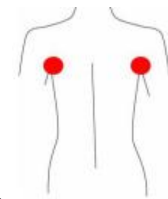


2. Circling by anticlockwise to stimulate the underarm lymph node.

3. On the back, by anticlockwise circling, can help to decompose the



fatness.

4. Using lymph drainage gesture, pull the toxin to the lymph node, can help improving the

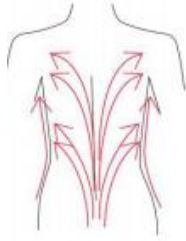

back curve.

ARM

Treated time: 30 minutes

BACK SIDE

1. From the inner elbow to the armpit, doing lymph drainage.
2. Circling to stimulate the lymph node.
3. By anticlockwise or come-and-go pushing, decomposing the fatness.
4. From the elbow to the outer, doing tightening gesture.



www.mychiway.com


Facade

1. From the elbow to the outer, doing tightening gesture.

Abdomen

Treated time: 20 to 30 minutes

1. From small to large, by clockwise direction, circling around the navel, can help peristalsis of the large intestine.

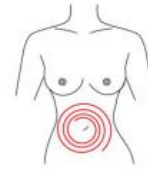

2. From small to large, by anticlockwise direction, circling push down.



3. From the belly and muscle texture, pull to the groin.



4. Through lymph drainage gesture, taking the toxin to the place of groin.



Attention:

1. Pregnant women or women during in menses.
2. Epileptic.
3. Patients with malignancy.
4. Patient whose wound after operation has not healed up.
5. Acute inflammation or epidemical patients.
6. Whom with heart diseases or with heart pacemaker.
7. Who with kidney (gall-stone) disease.
8. Who was embedded metal object or silica gel.
9. Who in menses, birth control period, emiction incontinence period, or accepting the belly operation.
10. Whose body always takes much inner hot.
11. Who has the genetic hypersensitivity.

Notice:

1. Be sure to use the special ultrasonic gel.
2. Avoid knocking the us head.
3. For avoiding burning the head, during operation, please prepare enough ultrasonic gel.
4. Don't stay in one place, avoid treating on the bone.
5. Don't use disinfectant product to disinfect the us head. We suggest the wet cotton or dry towel enough.
6. Check the machine power cord if it is connected well.
7. If the machine will not use for a long time, please switch off it, and power off.
8. Please take off all metal decoration from the operator and patient's body.
9. If continuous using 1 hour, please pause the machine about 10 minutes, then using again.
10. During the operation, be sure not to accept other treatment.

www.mychway.com