

4in1 RF +7 Colour Photon+SONIC+VIBRA

SR-267

Use Manual



CATALOG

Preface	3
I Theory	4
II Features	7
III Functions	7
IV Specification	7
V Package Listing	7
VI How To Use	8
VII Operation Steps	9
VIII Forbidden Matter	13
IX Daily Inspection And Maintenance	13

Preface

Dear User:

Thank you for purchasing !Be sure to read this instruction manual carefully before use so that you can make full and correct use of your massager.Please keep this manual available for future reference after you have read it through.

This product does not belong to medical,please don't use for disease treatment

I Theory



a. RF

Bipolar molecules are mainly composed of water and some proteins, they can get into deep skin through PSD from the probes, but it doesn't hurt the cuticle. When the temperature reaches 39 to 40 Celsius, the collagen will be shrunk by heat, and then accelerate its production. The efficacy can go deep into derma 2-2.5mm, increase the activity of collagen makes your skin firm and elastic.

b. LED(Light Emitting Diode)

For treatment of acne lesions

It works efficiently without damaging the skin by using an especially sensitive LED light developed exclusively for beauty equipment.

(1) Red Light (415nm-450nm)



Stimulate the production of collagen. Collagen is the essential protein used to repair damaged tissue and to replace old tissue. Best for fine lines, reduce pore size, swollen after surgery.

If you press the for a longer time, it will appear 7 Colour :

(1) Red Light (415nm-450nm)

Stimulate the production of collagen. Collagen is the essential protein used to repair damaged tissue and to replace old tissue. Best for fine lines, reduce pore size, swollen after surgery.



(2) Yellow Light (580nm-590nm)

Provide enough oxygen to skin, tighten skin and lift skin

(3) Blue Light (430nm-450 nm) :

Acne bacteria contains Porphyrins; The blue light is able to kill Porphyrins due to the high affinity with this wavelength. The blue light furthermore has a calming action resulting effective on hypersensitivity.



(4) Green Light: (520nm-530nm)

- a, Biological effects: it is a natural color, a neutral and stable effect between the blue and yellow feature;
- b, treatment effect: regulation of skin gland function, improve oily skin, blackheads, acne, etc.;
- c, Application: oily skin.

(5)Purple Light:(590nm-600nm)

detoxification ,Make the skin more healthy



(6)Moon Green Light:(630mm-640mm)

Make skin more bright, more elastic

(7)White light:(500mm-520mm)

Make the skin better import, better absorption



c. Sonic & Vibrations

Ultrasonic utilizes sonic vibrations; it works one million times per second (1 MHz).

The ultrasonic wave could stimulate skin tissue and improve permeability of the membrane to promoted collagen tissue, reduced wrinkle and invisible appeared fine line, also lift the skin result for younger and healthy look. Ultrasonic generate frictional heat which accelerate blood circulation, muscle relaxation, and decompose the cellulite to break down exceed fat.

II Features

- a. Time Master combined with 4 technologies with RF, 7 Colour LED, Sonic & Vibrations.
- b. Daily skin care for body & face therapy. Stop time flowing through your skin.
- c. Handy Time Master that is easy to carry around, treatment fulfill everywhere at anytime.
- d. Enjoy skin care at home with in-office experience.

III Functions

- a. Effect a radical cure blain blain, and curb the youth beans and black of the original bacteria, improve acne scars
- b. Remove wrinkles, skin lifting ,skin Rejuvenitation.
- c. Improve skin, promote the formation of collagen, make skin younger
- d. Remove the red spot

IV Specification

- a. Input Frequency: 50-60Hz
- b. Voltage: 100-240V
- c. Power Consumption: 10W
- d. Output Frequency: 90KHz

V Package Listing


- a. 1 pc Main machine
- b. 1 pc Power Adapter
- c. 1 pc package box

VI How To Use

a.Before use

Please connecting the power cord with the machine,and clean your skin

b.Process

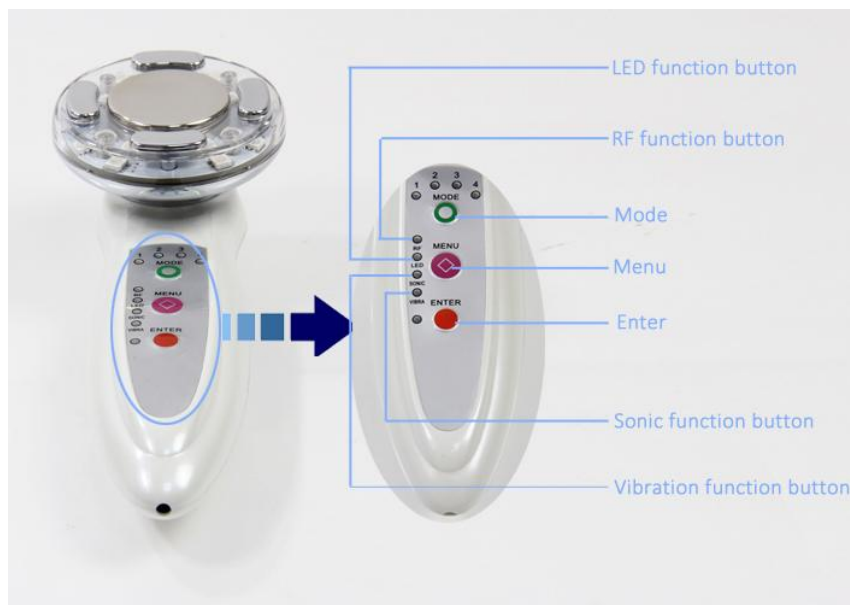

(1)RF

Press MENU for 3 seconds.Power on RF starts working.

Press MODE.Choose intensity from 1-4 accordingly.

After treatment,Press ENTER ,RF stop working.

(2)LED

Press MENU,choose LED.

Then press ENTER,LED starts working.

Press MODE ,choose Colours.

After treatment,Press ENTER again,LED stops.

(3)SONIC

Press MENU,choose SONIC.

Then press ENTER, SONIC starts working.

Press MODE ,choose intensity from 1-4 accprdingly.

After treatment,Press ENTER, SONIC stops.

(4)VIBRATION

Press MENU,choose vibration.

Then press ENTER, vibration starts working.

Press MODE ,choose intensity .

After treatment, Press ENTER, SONIC stops.

After all the treatment, Press MENU for 3 seconds, power off

The four function can work together.

VII Operating Steps


Once 2 to 3 days, 15 to 20 times a course of treatment

a. The leg(front):10 minutes


(1) From middle to two sides, it can active skin.



(2) From down to upper push, lift the skin.

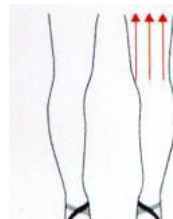


(3) From down to upper circling push.

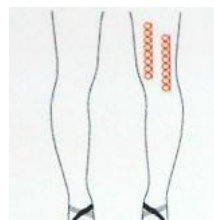


b. The leg(behind):10 minutes

(1) From down to upper push, lift the skin .

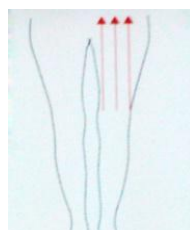


(2) From down to upper circling push.

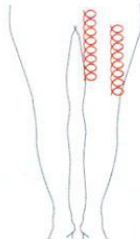


c. The waist: 10 minutes

(1) From down to upper push, lift the skin.

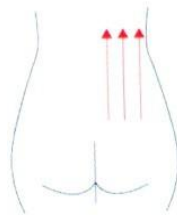

(2) From left to right push or from right to left .

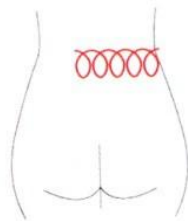

d. The arm: 5 minutes


(1) From down to upper push, lift the skin.



(2) From down to upper circling push, Stimulate the production of collagen.




(3) Outward from the arm looped, promote the effect of light, whitening, tender skin .

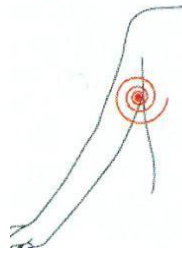


e.The hip: 5 minutes


(1) From down to upper circling push, active the skin.

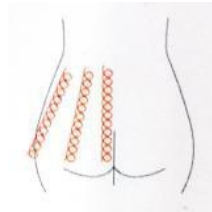

(2) From the crotch to the middle push, lift the hip, make Better hip shape .

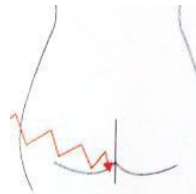


(3)From down to upper push, lift the hip stimulate the production of collagen.

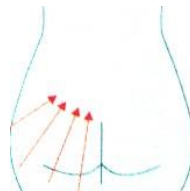


f.The waist:10 minutes


(1)push it as the line,active the skin on waist.

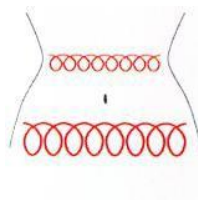

(2) In the middle of the waist, circle push around your belly button.

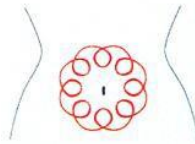

(3) From the waist of the lateral to the middle, Push it as pointed direction

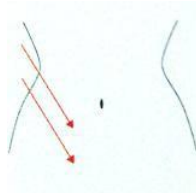

g. The face: 10 minutes.

Push it as the pointed direction.



h. The eye: 5 minute The lips: 5 minutes The neck: 5 minutes

Push it as the pointed direction.



VIII Forbidden Matter

Please noticed that this item is not applied to people as following:

- a. People with heart disease or high blood pressure, or configured cardiac pacemaker
- b. Patients with acute inflammation, asthma, deep vein thrombosis, theyroncus, cancer
- c. Pregnant women
- d. People with hemorrhagic disease, trauma or who is bleeding
- e. Medical plastic parts, or parts with artificial in filling inside
- f. People with metal inside the body (include metallic tooth)
- g. Patient with skin inflammation or with edema
- h. People with an abnormal immune system
- i. Numb or insensitive to heat

IX Daily Inspection And Maintenance

- a. If the head stick to a cosmetics or dirty, please use soft cloth gently wipe dry;
- b. Please don't tie line or winding, otherwise might well contact ;
- c. Thinner, volatile oil will be damaged, please don't use;
- d. Please do not placed in the sunlight, dust, more humid place (It is not waterproof);
- e. Please put it out of reach of children .

