

[Guarantee Certificate]

Product Model:
Purchase Date:
Purchase Place:
Sales Dept:
One Year repair Guarantee.

- The device will be repaired for free charge during the guarantee period. Please provide this Guarantee Certificate.
- The Guarantee condition not suitable following list:
 - Use the device not following this user manual.
 - Incorrect or unauthorized repair
 - Accidental damages during the transportation
 - Drop the device on the floor.
 - Damaged the device by fire, earthquake, weather such as flood, lightning and high wind.
- This Guarantee Certificate only suitable for China Mainland
- Please keep this Guarantee Certificate in safety place, will not duplicated again.

5. It is better to use the device after taking a bath.
6. The correct hold the device is keep upright on front side, insert the device as shown in the picture 9, five fingers grasped firmly.



Picture 9

Description of parts operated by the device



The device needs to move from the front of elbow to the underarm in "N" trace along the same direction with duration at 10 minutes.



The device will operate from the outer side of belly to the inner side in "N" trace along the same direction with duration at 10 minutes.



The device will operate from the knee to the upper side of leg in "N" trace along the same direction with duration at 10 minutes. The operation for front, rear, left and right of leg is the same.

Safety Notice

- Before use this device, please read "Safety Notice" in the manual and then use it correctly.
- Please take off your watch, necklace, bracelet and other metal decorations.
- Please use the device according to operating steps.
- Do not disassemble or transform the device by yourself, for repair, please consult the repair shop.
- Do not use this device beyond the description of this manual.
- Do not place the device under sunlight, near to brazer and bathroom or other place with high temperature and humidity.
- Do not place the device near to hot electric devices, which will cause fire, electric shock or other accidents.
- The power adapter can no be placed under wet condition, which will cause fire, electric shock or other accidents.
- Do not use the device once power wire and adapter has destroyed, which will cause fire.
- Keep the device far away from Children and Physical fault person.
- The device must be close to the body parts which need to be losing weight. Please stop to use it once you feel electronic stimulate too strong.
- The CV ultrasonic function must be match with slimming essential oil or lose weight cream together.

Specification and Package

Name: RF Radio Frequency Cavitation for Slimming
Input Voltage: AC100-240V
Output Voltage: DC12V 2A
Power: 24W
Net Weight: 222g
Product size: 114 x 90 x 100mm
Gift box size: 195 x 150 x 98mm



RF Radio Frequency Cavitation for Slimming User Manual

Please read this User Manual carefully before your use this device.

Product Introduction

RF Radio Frequency Cavitation for Slimming is popular in Japan; the main functions are Red LED light, RF Radio Frequency and CV ultrasonic. This is rechargeable type, can use it in wireless.


LED function: Red light with 625nm wavelength, this is high purity and strong light source, with uniform energy density, which has a significant effect in skin care and health treatment. Red light can enhance the vitality of cells, promote cells' metabolism, and make the skin release a huge amount of collagen and fibrous tissue for self-filling. Red light can suitable for various kinds of skin, promote blood circulation and regenerate the collagen to keep your skin more smooth, not only can improve the skin dryness, but also can improve hair loss and alopecia disease.

RF radio frequency function: Triangle Metal Conductive tips radiate radio wave of 500 KHz vibrations per second, this high frequency vibration will generate thermal energy to shrink collagen and attain the effect of tightening skin, improving wrinkles and reducing weight through continuous regeneration and recombination as the time goes by. Besides, the RF intensity of this device have 4 levels from weak to strong, so people can adjust this intensity according their own feeling.

CV ultrasonic function: The physical vibration with 330 KHz, make every body parts to vibrate and massage skin cells. Promoting the blood circulation, improving the nutrition of skin tissue, also keep the skin rosy and elastic. The ultrasonic can smash fat mass and make lymphatic detox. Thermal effects of ultrasound can change blood vessels and metabolism, strong circulation of blood, reduce cells excitement, so that analgesic effect. Ultrasound can enhance catalyze and accelerate the role of metabolism, which make the PH value change to alkaline, it is very good for drugs penetration, improve drug disinfection.



Product Parts

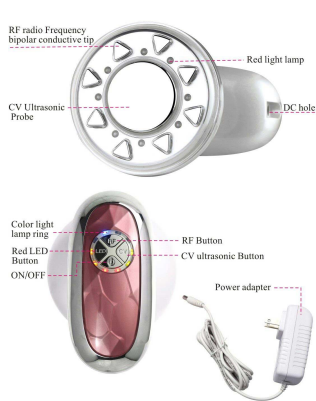

Usage

- Insert the Power adapter into Power socket.
- Put conductive gels or cream on the body parts.
- Press ON/OFF button, all color lights bright one by one, and then goes out.
- Press RF button, RF start working and warming, this is burning body fat function. You can adjust the RF intensity according your body feeling. Press RF in first time is first level, total 4 levels, from weak to strong, people can adjust levels according body parts, for example: there are too thick body fat at abdomen, so we can use 4th level which is the strongest thermal energy can burning body fat effectively; there are too small fat at arms, so we can use 1st level which is the weakest thermal energy, but it is enough for arms.
- Press CV button, this Ultrasonic function, can blasting body fat, ultrasonic can penetrate the skin and decompose the body fat.
- Press LED button, Red light with 625nm start working. Red LED can promote metabolism and remove wrinkles. Press LED twice the red light will flashing.
- Suggest use 10minutes for each body parts, the device build in auto off function is 10minutes.



Picture 5: CV Ultrasonic blasting body fat



Picture 6: Press LED buttons



Picture 7: Press LED button, LED light promote metabolism



Picture 1: Press ON/OFF button



Picture 2: Press RF button




Picture 3: RF thermal conductive probe



Picture 4: Press CV button

Special tips

- This device have three function: RF, CV and LED, you can use those three functions separately or in combination as following type: RF+LED; RF+CV; CV+LED; LED+CV+RF
- This device is applicable for the following body parts: Fore of arm, belly, upper leg, lower leg, back, waist, and hip (please see picture 8 shows with circle).
- This device is not suitable for following body parts: face, neck, clavicular parts, elbow, wrist, joint parts of body, knee, ankle, hindbrain, cervical vertebra and foot. (please see picture 8 shows with cross).
- Operating duration and times:
 - The size of each body part should within 30 x 20cm, move the device from up to down, left to right several times, around 10minutes.
 - You can choose three body parts every day, each body parts need 10 minutes for massage, so total is 30minutes. In same body part, you only can massage 2 times per week. It is better to massage it every 3 days.
 - The operation of this device needs to be matched with special slimming essential oil or lose weight gels for achieving an optimal effect.
 - After finished body slimming massage, please clean up the body, and apply with body lotion after 30 minutes.

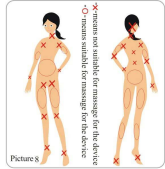

Picture 8: The size of each body part should within 30 x 20cm, move the device from up to down, left to right several times, around 10minutes. You can choose three body parts every day, each body parts need 10 minutes for massage, so total is 30minutes. In same body part, you only can massage 2 times per week. It is better to massage it every 3 days. The operation of this device needs to be matched with special slimming essential oil or lose weight gels for achieving an optimal effect.

WARNING
Please stop use the device immediately when you feel temperature too hot (especially in mode of RF level 4), avoid to burn your skin.