

Your profession family SPA care expert




COOLING AND WARMING PHOTON BEAUTY APPLIANCE


User Manual

Please read this user manual carefully and keep it properly before using the device

COMPONENTS

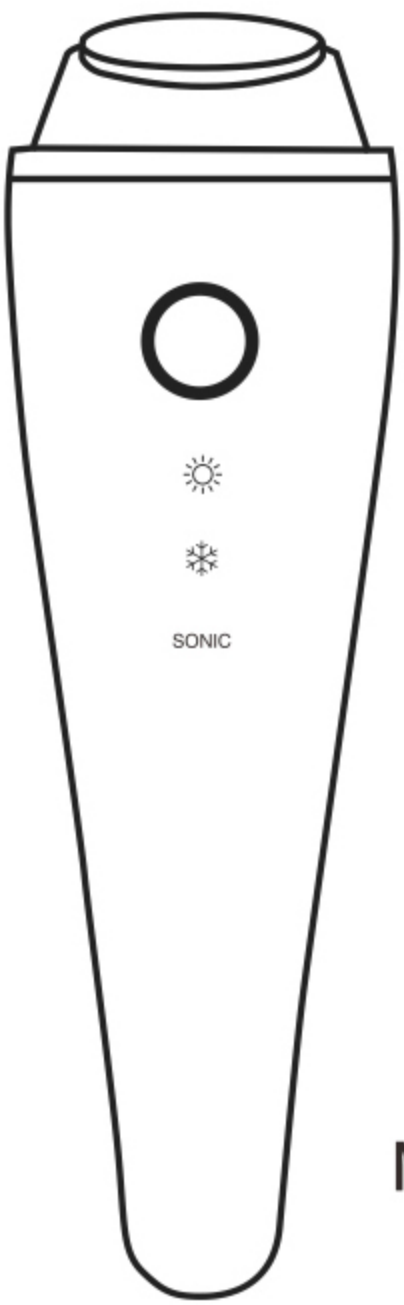

Main body

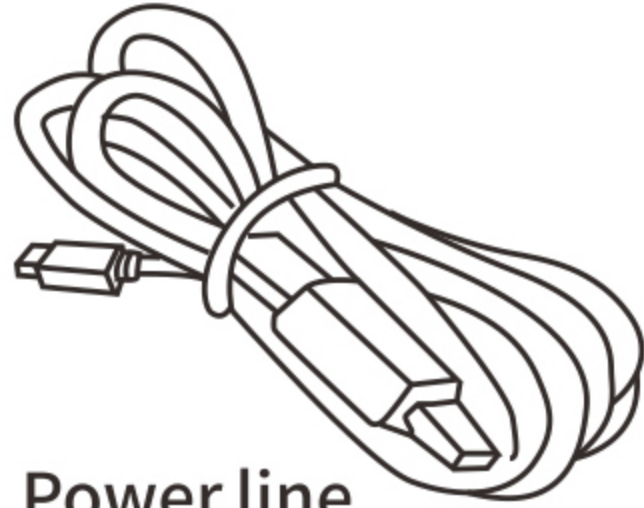

Power line

Please recharge the device before use

When charging: Device begins with red light, the light turn into green, charging stops if you put the cable off from the device.

Working hours

After fully charging, warming function of product's two terminals can use about 60 minutes, The larger one can use 40 minutes (the lithium battery capacity is 3.7v/900mAh)

Please make sure to charge before use, do not use it during charging, USB cable takes longer for full charging than power adapter.

Product specification

Product Model: BZ-0505

Weight:161G

Product measurement: 16.2*4.8*4.8CM


Input power: 5W

Charging time: about 2 hours


Battery life time: Warming 60min, cooling 40min

FUNCTION INTRODUCTION

Warming massage: Using 42°C constant temperature to dilate pore, promote the blood of circulation, accelerates metabolism, eliminates toxins.

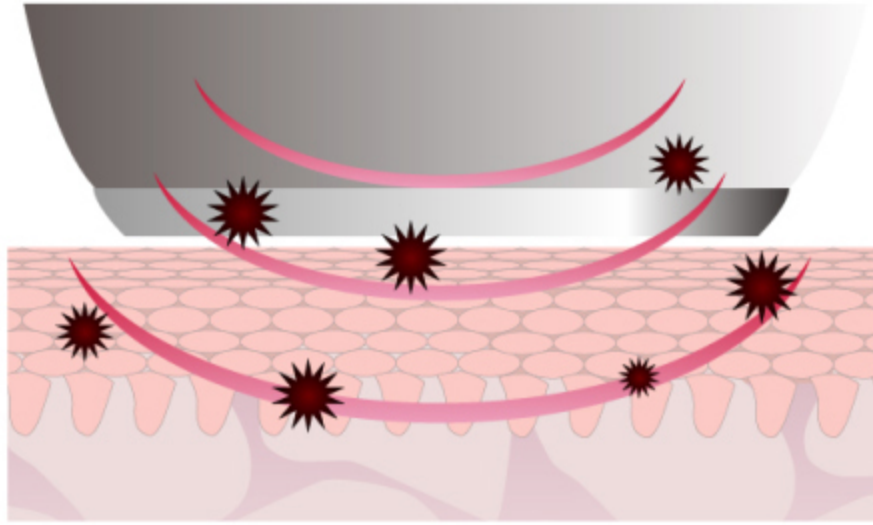

Cooling massage: Using 10°C constant temperature to shrink pore, lock-in water of skin, eliminates facial swelling when getting up in the morning.

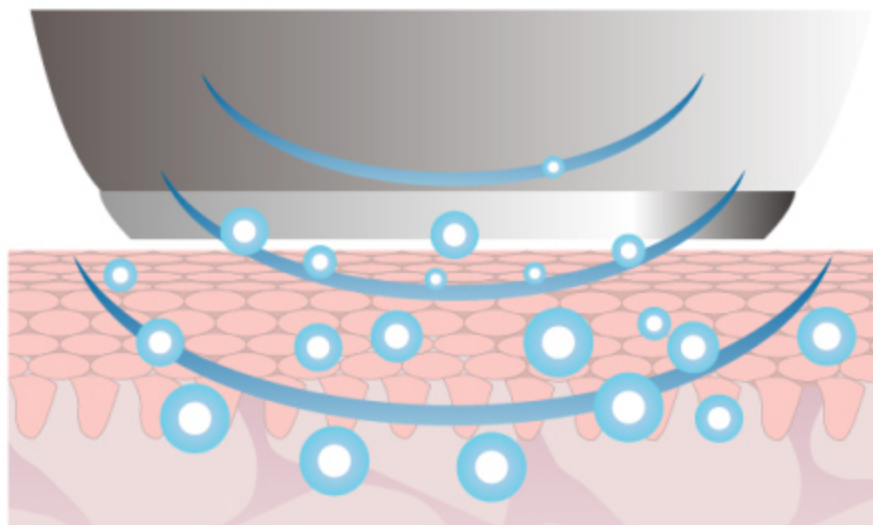

Vibration massage: Effectively stimulate the skin, the skin elastic.

LED red light care: Deep in the bottom of skin, improve blood circulation, promote skin metabolism, activates cell and collagen, eliminates toxin and beautifies skin, improves pigmentation, and recovers skin elasticity.

LED blue light care: Anti-inflammatory, reduces acne, relieves facial flush, good for acne treatment.

Safety Tips

- Please be sure to turn off when not using the device.
- Do not disassemble or damage randomly.
- Do not clean with chemicals such as alcohol, benzene after use.
- Do not place the device near hot or fire.
- Do not store the device where humidity is high or put it in water.
- Do not store near a microwave or high-pressure equipment.
- When charging is completed, disconnect the cable from the product.

CAUTIONS

- Please read this manual thoroughly before using the machine.
- Use after charging.(Do not use while charging)
- Ensure to use the provided adapter only.
- If you have a problem with the skin in use, please stop using it immediately.
- Do not use this device while bathing or hot shower.
- Avoid any place where humidity is high.
- Please be cautious of burn or chilblain.
- Do not drop the device or give it a strong shock.
- If damage is suspected or confirmed in the device, please do not use and contact your distributor or head office.
- Not recommended for pregnant women.
- Do not let the kids and the elderly use the device alone.

HOW TO USE

- ▲ Please check if the device is fully charged

Mode 1 (skin care)

Hot massage +LED red light+vibration

Please put good amount of cosmetics such as essence, facial cream, etc on your face before use.

Press the power button with duration 4 seconds and hot massage mode begin, promote the absorption of cosmetics. (It' s better to match facial mask)

Mode 2 (skin care)

cold massage + LED blue light + vibration

After cosmetics absorption , press the power button again, and the cold massage begin, shrink pore, make skin more elasticity.

Mode 3 (eye care)

Vibrating massage + warming massage

Finish facial care, press power button again, shift to eye care, relieve black rim of the eye and desalination eye bags effectively.


- ▲ The device will auto power off after working 10 minutes, if you want to use it again, please press power 4 seconds and the power will turn on.

After finished the device, please turn off the power.

The product can not be washed by water, but can use twist dry cloth to clean the remain dirt.


HOW TO MASSAGE

When using the device, please put the working probe on the skin lightly, then follow arrows to remove slowly.

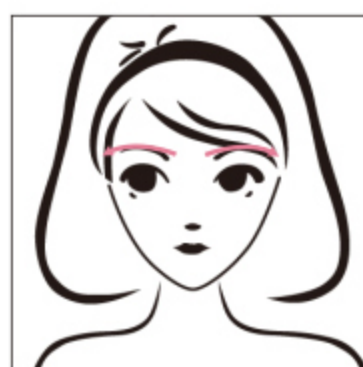


Forehead operation

Trace along the middle of eyebrows to the temple, from bottom to up, remove it gently.

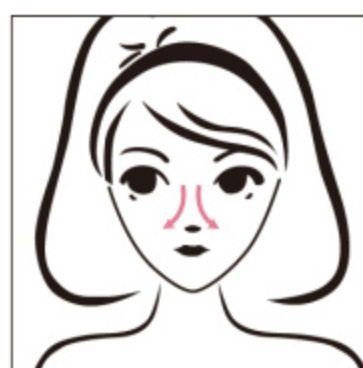


Nose bridge operation

Trace along the bridge of the nose, from up to down, remove it gently.

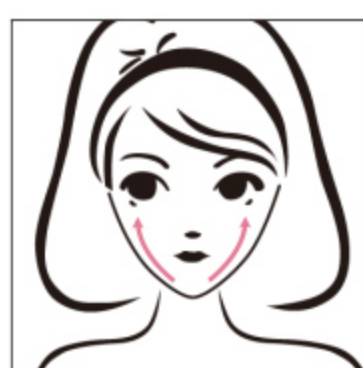


Cheek operation

Trace along cheek, from bottom to up and inner side to outer side, remove gently.

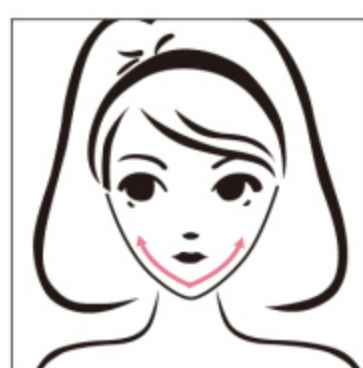

Chin operation

Trace along chin, from right to left, left to right, remove it gently.

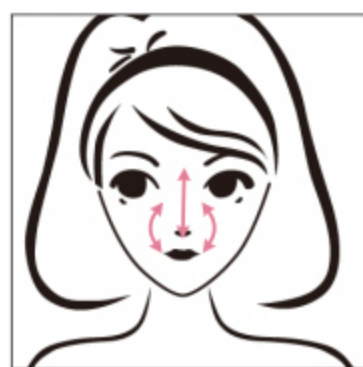

T part operation

Trace along wing of nose to nose bridge, from up to down for one to two minutes.



Lymphatic detox operation

Trace along side face, from down to up then go along the lower ear and neck, remove it gently.