


Operation skill

A.Face

Time: 10minutes

Put essence on the skin and then follow the skill as the indication arrow

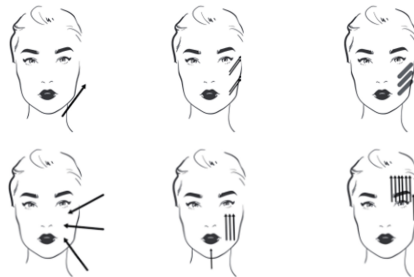

B Eye' s around

Time:2 minutes

Put essence on the skin and then follow the skill as the indication arrow



C Neck

Time: 5minutes

Put essence on the skin and then follow the skill as the indication arrow



D Arms

Time: each arm 5minutes

Put essence on the skin and then follow the skill as the indication arrow



E.Waist

Tiem: 10 minutes

Put essence on the skin and then follow the skill as the indication arrow



F.Butt

Time: 10minutes

Put essence on the skin and then follow the skill as the indication arrow



G Thigh and leg

Time: 20 minutes

Put essence on the skin and then follow the skill as the indication arrow



Time : 20 minutes



Time: 10 minutes



The patient who have cancer even past,especially skin cancer and disease pre-state  
 Suffering heart disease and other multiple severe disease  
 The person whose immunity brought down by AIDS and other immune function diseases,or is using immune to the preparations  
 The disease of recurrent herpes simplex and heated stimulated disease  
 The repair part had suffered lumps, psoriasis, eczema, rashes and other active skin disease  
 The repair part had surgery , laser resurfacing and the patient who needs enter the full course of treatment of deep peeling  
 Pregnancy and lactation  
 Allergic disease  
 If the user had accepted other beauty treatment and had suspected side effect,please consult the physician before using.

**Don't use on the following body locations:**

- 1、 Stomach when having menstrual
- 2、 The part after plastic surgery
- 3、 Wound
- 4、 Part of the acne and pimple inflammation
- 5、 Part that have dermatitis caused by cosmetics

- 6、 Part of swelling after basking and inflammation
- 7、 Part of pain and perceived disability
- 8、 Part which implanted metal , plastic , silicon and so on
- 9、 Part of abnormal spot caused by an itch , a fever , physical stimulation and so on

**Attention:**

This product is non-waterproof ,should not be flooded , splashed or contracted with any other fusible paste,otherwise it will cause an internal fault.  
 Using or store the product should not surround by hot or heating device,so as not to cause deformation or internal fault.  
 Don't put pressure and a strong impact on the product on it or make the contact deformation,otherwise it will cause the product's fault.  
 If the user has a discomfort half-way or under way,just stop immediately and go to a doctor.

**Operations:**

- Function:  
 A——RF heating  
 B——Needle-free Mesotherapy + LED light  
 H——the highest output intensity  
 M——mid-range output intensity  
 L——lowest output intensity


**Operation:**

Connect the power and the host is in standby,light A and B is flashing

Install the operating unit to the host,and the display window will show corresponding function code and flash.

Press the "Start" and indicator "L" flashes,press the "UP/DOWN" to adjust the output intensity.

Press the "Stop" to stop output.  
 Technical parameters:

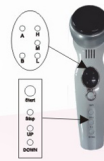

Size	18.3x5.9x7.68CM	Waterproof	Non-waterproof
Operation	remote control	Color	white
Input voltage	110-240V	Type of needle-free	high
Type of heating	Deep heating of RF	Wavelength of light	415-450nm
Operating voltage	DC12V 2A	Heat	45-65℃

(The above information is for reference only,if a minor change, please subject our product.)

**Disclaimer:**

- ※ Due to the complexity of body structure and physical difference,above effect is different from person to person.
- ※ Deep heating of RF can accelerate blood circulation ,but cannot replace the drug therapy.Prohibit stopping the drug without doctor's permission.



**Hand-held RF Skin Rejuvenation Device**

**Hand-held RF Skin Rejuvenation Device**

**Description:**

Hand-held RF Skin Rejuvenation Device adopts the deep heating of the RF, needle-free mesotherapy and light technology, it makes them a fantastic combination.It devotes to make the subepidermal cells movement active,enhance the effect of skin management, improve skin's elasticity, astringe pores , repair wrinkle and improve the permeability of cosmetics,so that it can help us rejuvenate the skin.

**Application of Principle:**

Deep heating of RF technology:  
 The deep heating of RF effects the electronics of tissue with the biological response of polarization electronic mobility by way of human tissues,resulting the electronics are formed as the molecules twist and crush against each other so that bioenergy is produced,thenby bringing skin to warm in deep to stimulate the collagen produce contraction immediately,to stimulate the secretion of new collagen to fill the gap of losing collagen atrophy ,and rearrange to rebuild skin soft frame ,and ultimately firm up skin ,remove wrinkle,restores the skin's elasticity and luster.

- Facial contours of the sculpture,such as face-lift , reduce a double chin ,diminish fine lines and laugh lines, and tighten loose skin, instead of and beyond the surgical face lift.
- Evacuation of lymph, eliminate excess water and toxins, remove eye pouch and dark circles, astringe pores and acne.
- Dispel false authenticity wrinkles and repair stretch marks.
- Soften the horniness ,reduce melanin and purify and whiten skin .
- Systemic anti- aging,tighten the modification of neck and back ,shape the waist , the abdomen's surrounding and the legs.

**Needle-free mesotherapy:**

Needle-free Mesotherapy is a combination of mesotherapy and electroporation.In Europe which has advanced cosmetology the mesotherapy without injection is very popular.Its feature is the infiltration is strong,the infiltration is from skin's surface to 1mm which is 10 times to ion's infiltration and to 4mm is about 27 times.It can really infiltrate the effective constituent into the skin,and it could play a high-performance,and thus the actual experience of the effect of peeling.

**LED light:**

Hand-held RF Skin Rejuvenation Device issued by the soft LED light which won't cause hurt to skin and also fully effective. The most suitable absorption wavelength is 415nm-450nm.LED light which with the calm skin effect can prevent and improve acne .

**The instruction of hand-held RF Skin Rejuvenation Device**

Sincere thank you for buying hand-held RF Skin Rejuvenation Device.The instruction is made for using this product safely .This product is very easy to use ,but according to the different method ,sometimes not only not get the full effect, but also dangerous.So be sure to read this instruction and understand fully before using this product.

**Danger:**The following instruction contains important safety information,which may lead to death or critical injury if not followed.

**Warning:**Injure may occur when the following instruction is not followed.

**Attention:**Injure and product failure may occur when the

following instruction is not followed.

**Danger:**Don't use the following medical treatment at the same time and prohibit wearing jewelry of precious metals,such as watches and rings.

- 1、 Body has medical planted with electronic instrument
- 2、 A heart-lung and life with electronic instrument
- 3、 Heart motor installed with electronic instrument of medical
- 4、 User who have been implanted with metal

**Warning:**

1、 Please don't do any repair , disassemble and transformation on the product or other related components,otherwise it will cause adverse phenomenon,the fire, burns and so on.

2、 If the product occurs an abnormal sound or abnormal phenomenon ,please stop using it immediately and repair. Be sure to use the power that is equipped,or it will cause an accident or malfunction.

Don't use when the user have the following symptoms :  
 The patient who are using a pacemaker or the implantable rectifier de-fibrillation  
 The repair parts had been implanted metal