


Mini Fractional RF Skin Tightening Beauty Device Model SR-W705 User Manual

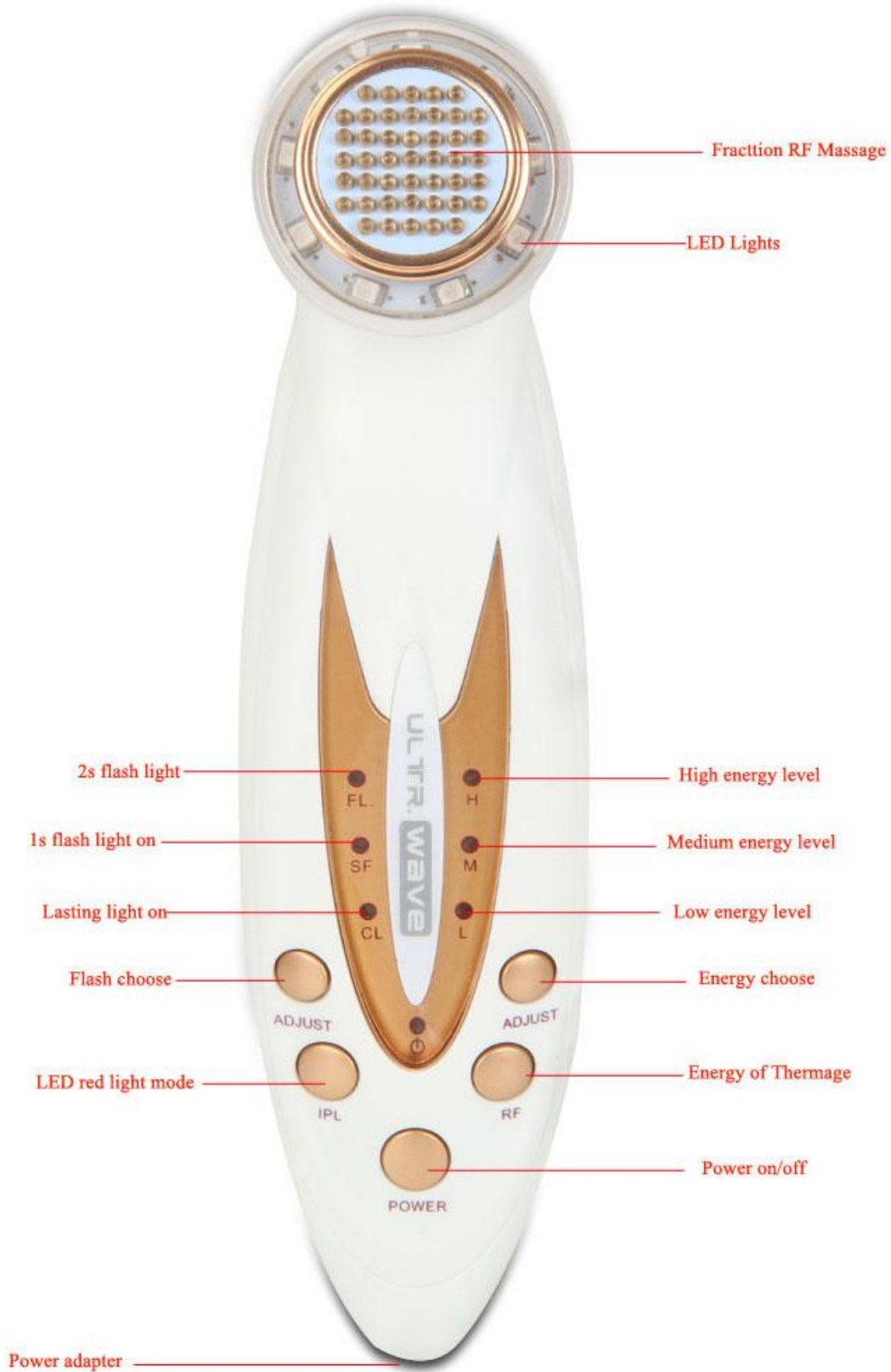

Using method (steps)

- 1, keep the maintenance area clean and tidy.
- 2, take appropriate amount of gel or extract liquid onto the face, neck and other maintenance areas.
- 3, press the POWER key then turn on machine.
- 4, press LED red mode selection button (IPL), then red light long-lighting mode starts. Then press ADJUST key adjusting requested modes.
- 5, press the RF key the low grade mode works. Then ADJUST key for required energy grade. (Before the first usage, please test the underside of your arms first, selecting suitable grades based on your acceptable degree of stimulation and temperature.
- 6, Different maintenance part with different and appropriate method.
- 7, after the maintenance press ON/OFF key then turn off the machine.

Attention!!

There will be skin fever or slight pain during the maintenance process; Some would feel skin rubefaction but will disappear in minutes. Please stop the maintenance if you suffer abnormal fever.

九 , commonly used maintenance ways

- 1, dredging lymph and back to young appearance.

Put the head of the beauty machine close to the proximal clavicle artery, from the inside to outside circling upward to the back of the ear; then pull down forcefully along the artery to the start point (proximal clavicle artery), repeat the whole movement for several times.

- 2, Remove double chin and eliminate

Put the head of the beauty machine close to the double chin, pull backward forcefully to back of the ear, then pull downward from back of ear to the proximal clavicle artery; and next the opposite direction, pull upward slightly from clavicle to cheekbone, then pull forward to double chin, repeat the whole movement for several times.

- 3, V face-shaping

a, Put the head of the beauty machine close to chin, circling forcefully upward to tragus side, and backward from tragus side to the back of ear, then pull downward from back of ear to clavicle; and next the opposite direction, pull upward slightly from clavicle to cheekbone, then pull forward to the chin, repeat the whole movement for several times.

b, put the head of the machine close to the corner of mouth, circling outward forcefully to ear side, along the ear side pull backward to back of ear, then pull downward to clavicle; and next the opposite direction, pull slightly upward from

clavicle to cheekbone, next pull forward to the corner of mouth, repeat the whole movement for several times.

C, put the head of the machine close to malar bone, circling forcefully outward to ear side, along the ear side pull backward to the back of ear, then pull downward from back of ear to cheekbone; and next the opposite direction, pull slightly upward from clavicle to cheekbone, then pull forward to malar bone, repeat the whole movement for several times.

4, Remove wrinkle near mouth and nose back to youth

Put the head of machine close to chin, push upward forcefully to philtrum based on the method of drawing letter C; then pull slightly downward from philtrum to chin with same method of drawing C, repeat the whole movement for several times.

Put the head of machine close to the wrinkles alongside the mouth, circling outward forcefully to the wrinkles (alongside the wing of nose), then pull downward slightly back to the wrinkles alongside the mouth, repeat the whole movement for several times.

5. Lighten nose shrink pores--promote collagen regeneration and increase elasticity

Put the head of machine close to the wing of nose, circling outward forcefully up to the bridge of nose, repeat the whole movement for several times.

6, remove eye pouch, eye wrinkle and black rim of eye

Put the head of machine close to eye pouch, at a distance about 1cm from the eyes pull upward forcefully to the temple, and repeat the whole movement for several times.

Rapid Repair Technique

Rapid repair awaiting repaired wrinkles, scars, freckles and other important parts

Put the head of machine close to repaired part for seconds (kind remind: 10 seconds with low-grade mode, 5 seconds with mid-grade mode and 1 second with high-grade mode; once stay for longer seconds your skin could be burned), then move to the end repaired part repeat the movements for several times. The procedure just likes seal on the skin, it ends up finishing all the repaired parts.

Specially attention:

Avoid from your eyeball. With low-grade mode to the tender skins around the eyes.

Deep leading-in techniques

Make deep leading-ins with facial treatment essence or nursing products containing with “active gene repair function”.

Put the head of machine close to the skin wiped with facial treatment essence or containing with “active gene repair function” nursing products, circling back and forth alongside the skin texture, make the machine head super close to skin during the movement with moderate i