

【6 in 1 body slimming machine】

6 in 1 body slimming machine



Operational Manual
[6 in 1 body slimming machine]

- ◆ Thanks for your purchasing our beauty device
- ◆ Please read this manual carefully before using, and keep this book in a safe place for future reference
- ◆ This product is designed for take care of skin and make skin beauty, and other usage or sales are unavailable

Contents:

Introduction of body slimming device -----	1
Safety and Warning -----	2
Accessory and parts -----	3
Operation method -----	4
Ultrasonic massage operation -----	4
Galvanic ion operation -----	7
LED therapy operation -----	9
How to use EMS -----	10
Specification -----	12

Introduction of body slimming device

This Device has six Functions:

- 一、 **Ultrasonic Function:** Ultrasonic Vibration with 1000000times/s, this vibration is very suitable for taking care body. There are two main effects on skin:
 1. **Mechanical Massage:** 1MHz Ultrasonic vibration with 1000000times/s, this high frequency and fine vibration, only hydrone can feel it, if you drop some water on the ultrasonic probe, you can see water spring and atomizing quickly. Sometimes our skin can not feel it, this is normal.
 2. **Warm Treatment :** Make skin temperature rise 0.5-1℃, so that nutrients can be quickly absorbed, accelerate metabolism, enhance tissue regeneration process, repair poor lymphatic and various capillaries, so make weight lost and slimming body, so ultrasonic can help lead in nutrient for skin effectively.
- 二、 **Positive Ion leads in function:** can lead much more nutrition deep into skin, make the skin cells fully absorb nutrient. Take care skin effectively.
- 三、 **Negative Ion leads out function:** match with clean cream, can deep cleaning skin, absorb residual cosmetics metal ions. Those residual cosmetics metal ions will become harmful sediment into skin; make the skin appear dark spot and rough aging.
- 四、 **Red light therapy:** 620-630nm Red light can anti-aging and beauty skin tighten skin pores and dilute dark spot, at the same time can provide energy of cell activation, promote metabolism and anti-aging. Promote blood circulation and proliferation of collagen.

五、 **Blue light:**455–470nm, Prevent bacteria growth of acne wounds, avoids skin inflammation, also can avoid leaving scars of acne. Special blue light can make skin shine and bright.

六、 **EMS Stimulation Function:** Using micro-current massage, under this mode, the current stimulated muscles, to make muscle exercise, and consume surrounding fat, sugars and carbohydrates. Under normal movement, muscle only can subsume 60-70 joules per second, but under this high-frequency movement by EMS, there will consume up to 3000-4000 joule per second, so that make your body sliming quickly and effectively.

Safety and Warning


Storage attention:

1. Do not use this device when you found power wire and power adapter is hot, damage or power adapter become flexible.
2. keep the device away from water or moisture place.
3. When you don't use the device, please cut off power.
4. keep away from children.
5. Patients, Pregnant women, and cardiac patient can not use this device.

Do not use this device for special people or some body parts, see following notice:

1. Eye ball, body parts which just use hair removal products/machines
2. The skin has inflammation due to strong sunlight
3. Dermatitis due to cosmetics
4. With pimple or inflammation parts, do not use galvanic ion and ultrasonic function. Only can use light function
5. The skin feel s itching or fevering
6. Some spots which caused by physical stimulation
7. Reshaping body parts
8. Heart disease patient (especially with pacemaker) an pregnant women
9. Other parts which is not including in operation manual

Accessory and parts

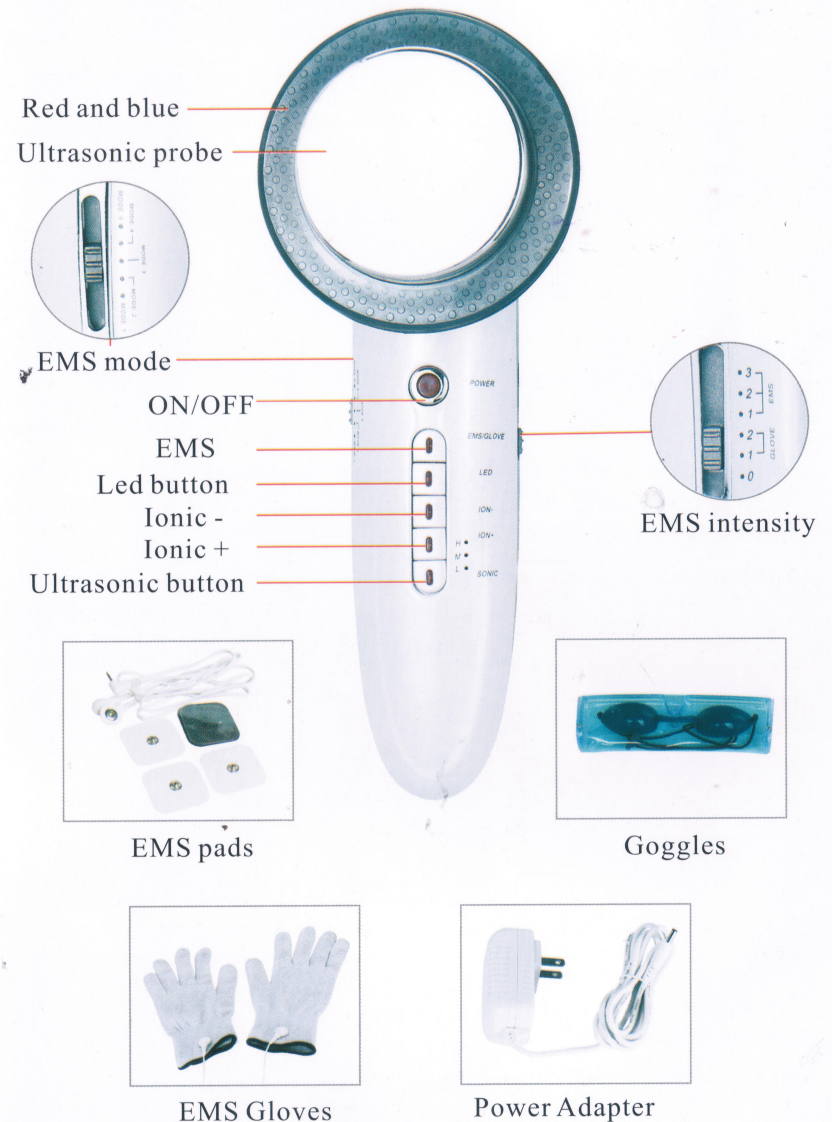


Application method

— Ultrasonic massage operation

(Before your using ultrasonic, it is very important to put some skin care cream or liquid on your skin)

- 1) insert the power adapter into socket
- 2) Press Power button on Device (Red indicator)

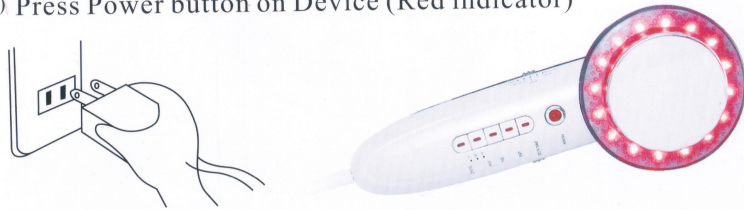

- 3) Press Sonic button, Sonic function start work, and intensity is Low, Press Sonic button secondly, intensity is High (please use Low intensity for thin skin)

- 4) Please put skin care cream/liquid or slimming cream/liquid on your skin, make you skin fully moisture,



- 5) Put some skin car cream/liquid or slimming cream/liquid on the probe of device.
- 6) your body according following method

Warning: Do not touch your throat, eye ball, parts of near eyes, Elbow, Knee and other joint parts.

4

Body Treatment

Following Body Treatment methods are only basic using operation, if you need more professional treatment, please ask beauty consultants

1、Massage Arms



Make arms drooping, relax muscle, and massage arms from wrist to the armpit.


Effective: shape arms and make arms slimming



Massage arms from wrist to the armpit by circle massage method.


Effective: reduce fat for arms

2、Massage waist and abdomen

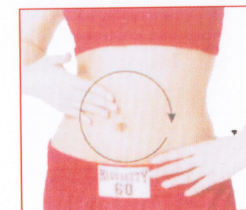

Massage abdomen according clockwise direction

Effective: tighten body skin, and remove body cellulite.

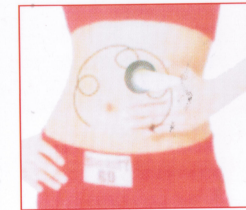

Massage abdomen according clockwise direction by circling method.

Effective: remove body fat and body wrinkle.



Massage waist according pictures show.


Effective: Tighten waist skin.

5


3、Massage hip (please keep the massage head contact with your skin directly)



Starting from the center of buttocks, massage your buttocks in a circling movement toward the back
Effective: reduce the cellulite of buttocks

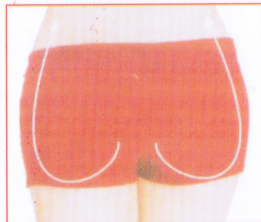

Massage buttocks from the end of coccyx to waist
Effective: Lift and shape buttocks

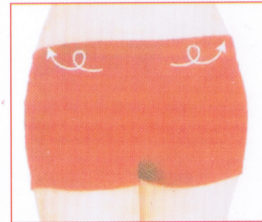

Massage buttocks from upper buttocks, apply circling movement
Effective: tighten waist and buttocks

4、Massage foot and legs



Massage from the low leg to thigh
Effective: lift leg skin and shape leg



Massage leg by circling method
Effective Remove the leg fat and wrinkles




Massage thigh from inner to outside
Effective: Tighten skin, prevent skin loosening

二、Galvanic ion operation


Operation for Lead out function

1. Use wet towel to clean the body.
2. Press power button (Red indicator).
3. Press +Ion button, you can hear "tick" sound, button show red indicator, that is Positive ion, start lead out function, it can deep cleaning skin.
4. Please hold the handle of device when you use ion function.

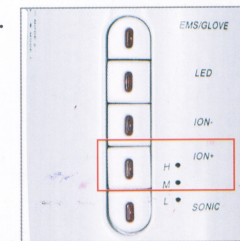

Operation for Lead in function

1. Use wet towel to clean the body.
2. Press power button (Red indicator).
3. Press Ion button by secondly, you can hear "tick" sound, button show green indicator, that is negative ion, start lead in function, it can import nutrition deep into skin.
4. Put some skin care gel/liquid or body slimming gel/liquid on your skin, make your skin fully moisture. (gel type skin care products will be more better)
5. Also put some skin care gel/body slimming gel on the sonic probe. Massage body and make skin care gels deep into skin.

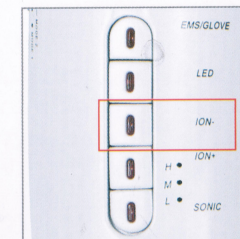


Note: The device with timer will auto off after 10 minutes, so if you need massage again, please power on it again.

We suggest you use this device 10minutes/time for each function if you use this device first time. Generally you can add time to 30minutes/time for each function.


You can use three functions together, total two modes:

1. Sonic + Red light or Blue light + Positive ion.
2. Sonic + Red light or Blue light + Negative ion.


Facial basic care

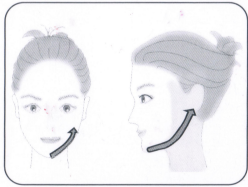


1. from chin to the back of ear

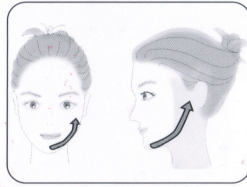


2. from chin to the front of ear

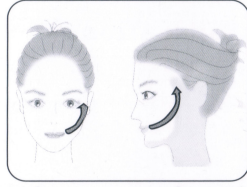


3. from corner of mouth to the upper of ear

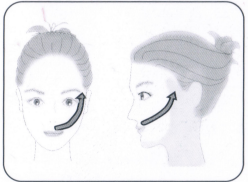

4. from bottom of nose to the near of eyebrow

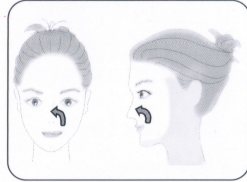

5. from wing of nose to bridge of nose

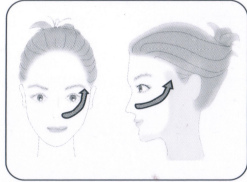


6. from nose to temple or forehead (Note: please don't stop in temple parts)

Facial Special care

Because of earth gravity, age increasing and fat accumulating, will make your skin loosening and pendulous, so appear double chins (Use this massage 10-20 times per stage)

Please give slight strength to ultrasonic probe when you do massage, lift your skin up according following arrow indication.

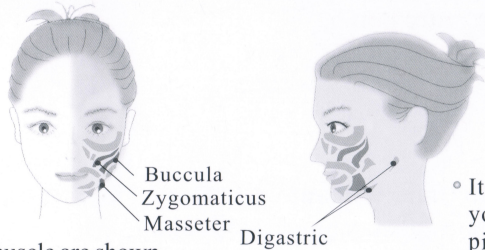


facial muscle are shown in the above drawing

It will be more effective if you massage face according picture showing

Operation to tighten face

Operate from the corner of mouth towards the area near eyes.

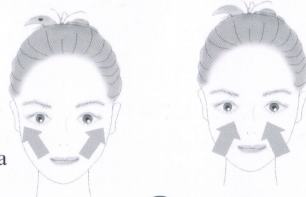

Operation to pull up the end of mouth

Operate from the corner of mouth in upward motion.

8

三、LED therapy operation

1. Press Power button, then press LED button. Blue light start work


Blue light remove acne.

a) Biological efficacy
High purity blue light can kill propionibacterium of acne, quickly healing the acne.

b) Therapeutic effect:

Shrink vascular, enhanced protein fibrous tissue, tighten pore, suitable for aging skin, and wrinkle skin and edema skin. Blue light have calming and resist acne function. Can adjust skin sebum secretion, remove acne.

c) Application: sensitivity skin and oil skin.

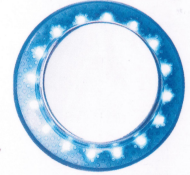

2. Press LED button secondly, Red light start work

Red light removes wrinkles


a) Biological efficacy.

Stimulated cells increased vigor, enhance blood circulation.

b) Therapeutic effect:

Lymphatic Drainage, Soften tissue, apply to the eyelid fat, eye bags, dark circles, remove spot, make skin moisture, elimination yellow face, remove wrinkle, accelerated circulation, promote skin tissue activation. Red light also can tighten skin.

c) Application: any skin.

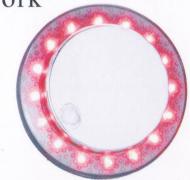

Note: can use three functions together. Two modes:

- I. Sonic + Red light or Blue light + Positive ion.
- II. Sonic + Red light or Blue light + Negative ion.

After using

1. Please Power off all buttons of device and then cut off the electronic power.
2. Use wet cloth to wipe dry probe, cleaning the probe.
3. please don't use some chemical liquid to clean the device.
4. Put device in box and keep it storage.

9

四· How to use EMS


1. Use EMS Pads

There are 5 modes for EMS function; people can choose any one according your feeling.

Suggest sticking Pads on following position;


BELLY TIGHTENING

1. SHUIDAO (20)
2. SHUIFEN (6)
3. GUANYUAN (18)
4. ZHISHI (53)

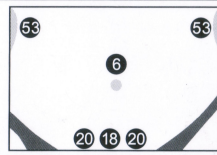


HIP TIGHTENING

1. HUANTIAO (54)
2. CHENGFU (62)

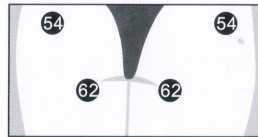

WAIST BEAUTY

1. TIANSHU (10)
2. TAIYI (7)
3. DAIMAI (17)
4. SHENYU (47)

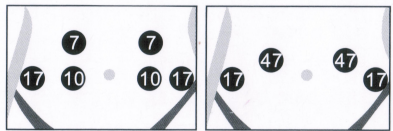


BELLY BEAUTY

1. GUANYUAN (18)
2. DACHANGYU (48)
3. ZHONGYUAN (8)




SKIN BEAUTY

1. FEIYU (44)
2. GANGYU (45)
3. ROUMEN (5)

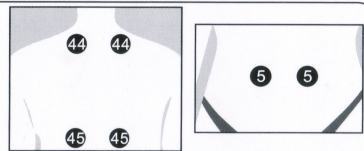


WEIGHT REDUCTION

1. DAZHUI (42)
2. YANGGUAN (52)
3. ZHONGWAN (8)
4. GUANYUAN (18)
5. JULIAO (22)
6. SHENYU (47)
7. CHENGFU (62)

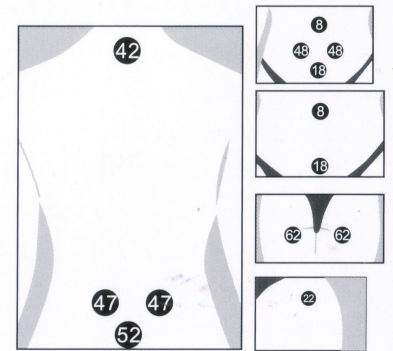

2. Use EMS Gloves

Use gloves not only can massage whole body, also can massage breast, stimulate blood circulation.

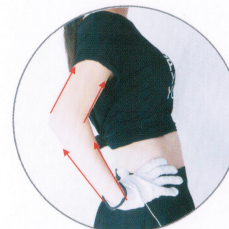


arms




face



abdomen

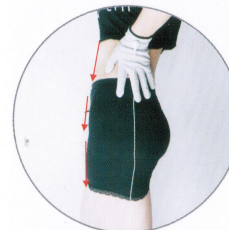

buttocks

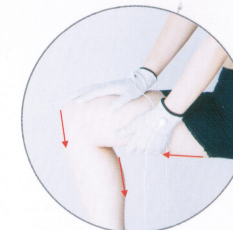

thigh



shank

Specification:

Name: 6 in 1 body slimming machinesliming device

Power adapter input: AC100-240V

Power adapter output: DC15V 800mA


Power: 6W

Frequency of Sonic: 1MHz

Red light wave length: 620-630nm

Blue light wave length: 455-470nm

N.W:191.4g

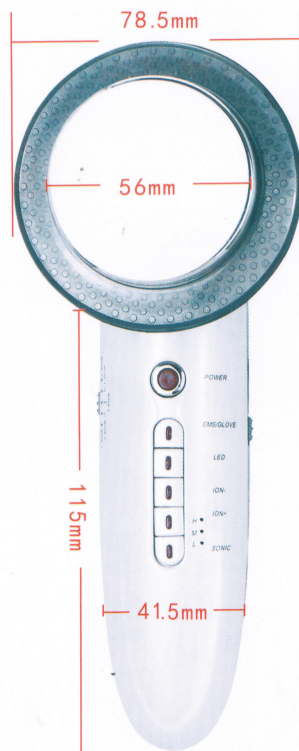

12

[Guarantee Certificate]

Product Model:

Purchase Date:

Purchase Place:

Sale Dept:

One year repair guarantee

1. Please use the device according user's manual, repair the device for free during the Guarantee period (please return damaged device and this guarantee certificate)

2. Repair fees will be charged if the device damaged due to following reasons:

- (1). Operate the device not according our user's manual
- (2). Tear down and repaired by customer self
- (3). Destroyed during the transportation or drop device on floor.

(4). Damaged due to fire, earthquake, bad weather such as flood, lightning an high wind.

3. This Guarantee condition only suitable for China Mainland.