

USER MANUAL



MS-12C2

Contents

1. Brief Introduction	2
2. Advantages	2
3. Working Principle	3
4. Effects	3
5. Indications	4
6. Contraindications	4
7. Matters Needing Attention After Treatment	4
8. Detailed Operations	5
9. Technical Parameters	5
10. Packing List	6
11. Precautions for Use	7
12.FAQs	8
13.Operational Diagrams	9

Brief Introduction

Thanks for choosing our beauty equipment for fat blasting and body slimming. This equipment applies 60K ultrasonic cavitation technology. It can effectively reduce local flab and obesity. By adopting high-tech, it manages body in a safe, comfortable and effective way. Besides, it can be operated easily, conveniently and swiftly, and solves local obesity problem for people pursuing beauty, and can replace weight reducing and body slimming programs in beauty salon. It's such an amazing device that you can reduce weight easily at home. Moreover, this equipment works safely and effectively, and requires no injection, no medication and no operation. It has no side effects, and is operated externally throughout, and makes people feel comfortable. MS-12C2 is a professional body slimming equipment, which needs to be operated by people with professional training. All the people are required to read this manual thoroughly and follow the instructions specified hereby strictly before operation.

Advantages

- 1.Ultrasonic cavitation strongly blasts and removes fat, and has replaced and surpassed liposuction, which brings good news to obese people.
- 2.It adopts the most advanced cavitation technology in the world.
- 3.It's comfortable, painless, noninvasive and injury-free during treatment, and there is no convalescence, and without affecting people's normal work and life after treatment.
- 4.No consumption, low cost, but with quick return.
- 5.Unevenness, bleeding and swelling will not occur.
- 6.It's compact, easy to carry, to learn and to use, and can be both used at home and in beauty salon.
- 7.Breakthrough: It combines together heating, vibration and ultrasound. Compared to traditional ultrasound fat-dissolving, this equipment has a more obvious and evident effect.

8.It heats up deep skin, vibrates strongly, accelerates blood circulation and rapidly discharges fat from body.

9.Ultrasound aims at body parts having difficulty in fat reduction, and will bring amazing effect.

Working Principle

It acts on body by utilizing a strong sound wave head of 60,000HZ, which goes deep into fat layer vibrating rapidly deep and obstinate fat. High frequency vibration can cause friction among fat cells, and effectively consume heat and cellular water, thus contracting fat cells. At the same time, one-sided vibration causes strong impact to fat cells, which generates cavitation effect. Small bubbles formed by cavitation vibrate constantly, grow or burst instantly along with the vibration of surrounding medium. When hollows that grow inside and outside of those fat cells get burst due to pressure changes, fat cells also burst, thus dissolving triglyceride into triglyceride and free fat acid which achieves the ideal effect of fat removing. 85% of fat cells decompose and discharge from lymphatic system, while the remaining 15% are translated into carbohydrate via liver which will be consumed in daily activities. Breakthrough: It combines together heating, vibration and ultrasound. Compared to traditional ultrasound fat-dissolving, this equipment has a more obvious and evident effect.

Effects

- 1.Strongly blast fat, rapidly eliminate fat, dissolve fat and sculpt beautiful body curve.
- 2.Accelerate detox and metabolism, and increase skin's elasticity.
- 3.Blast fat, remove flab, improve local fat accumulation, swim ring-like belly, thick leg, etc.

Indications

1. People who always sit without exercise, or with jiggly arms, swim-ring like belly, thick leg, or flab on waist and abdomen.
2. People whose lower part of body is obese, or with pear-shaped body shape.
3. People with local obesity, fat accumulation or who looks ordinary in clothes.
4. People with unsymmetrical body shape, or who are overweight.

Contraindications

1. People in pregnancy, lactation, menstruation, or surgical recovery.
2. People with serious diseases, such as hypertension, heart disease, diabetes, severe thyroid, malignant tumors, etc.
3. People with skin disease, infectious disease, or whose skin is in sensitive period.
4. People with wound, or who just had liposuction.
5. People with severe varicose veins or tumor.
6. People in allergic period, or who have severely sensitive skin.
7. People who are overly aging.

Matters Needing Attention After Treatment

1. Keep warm, and don't eat cold food, and avoid exposing to windy environment or catching a cold.
2. Take shower after 4 to 6 hours.
3. Drink more warm water, and replenish moisture to speed up metabolism.
4. Refuse eating and drinking too much, as well as staying up late.
5. Avoid enjoying sauna, hot springs, or doing strenuous exercise within 7 days.

Detailed Operations



- Remark: It's advised to set the energy level between 3 and 7. For first time user, it's suggested to start with the lowest energy level and increase it gradually.

Technical Parameters

Size:	230(L)*190(W)*95(H)mm
Package Size:	310(L)*260(W)*102(H)mm
Power Supply Input:	110-240V
Power Output:	Max. 30W
Frequency of Cavitation:	60KHz

Packing List

Main Machine ×1

Handle x1

Power Cord ×1




Main Machine



Handle



Power Cord

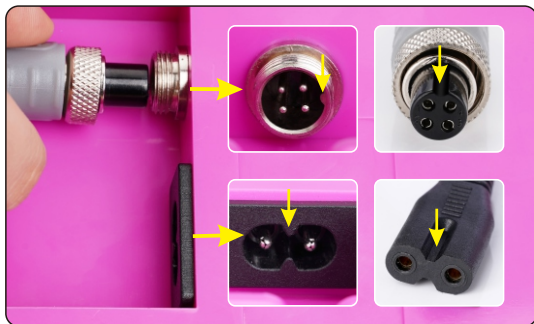

Connect handles to the corresponding sockets

Precautions for Use

- 1.Clean the operating handle with normal saline after use and keep it properly so as to extend its service life.
- 2.Plug with ground pin must be used, and power socket that already factually grounded must be ensured before using the equipment.
- 3.To make sure equipment's voltage is adaptive. If voltage of local power supply is unstable, we suggest users add voltage regulator with matching power between the mains supply and the equipment.
- 4.To guarantee curative effect and normal service life of the equipment, please uniformly use specified parts provided or suggested by original manufacturer.
- 5.Don't place the equipment near a strong heat source since this may affect its service life and normal use.
- 6.Please remove all the metal objects from body before treatment so as to avoid unexpected situations which may affect curative effect.
- 7.Please restrain from using the equipment aiming at eyes, thyroid, parathyroid, testicles, pregnant woman's abdomen, pacemaker, etc.
- 8.People who are suffering from illness should use with caution unless he gets permission from a doctor.
- 9.Please turn off power switch of the equipment if no one uses it, and ensure main power is off after someone used it and before he left, thus safety of the equipment can be guaranteed.
- 10.It's advised to stop taking other weight-loss medicines for 1 to 2 months before treatment if you are taking them now. The course of treatment shall be extended if you want to have it started immediately.
- 11.Don't use the equipment with an empty stomach, and wait for at least 1 hour after a full meal before starting the course for weight reduction.
- 12.The equipment should contact skin fully so as to avoid uneven heating when in use.
- 13.Start with the lowest energy level when just started, and increase the level gradually after the client got used to it.
- 14.When using this equipment, parts to be treated must be kept moist, and dry skin treatment should be avoided.

15. Avoid contacting skeletons, heart and thin-skinned parts as much as possible when in use.

16. Using the equipment and training the operators in strictly accordance with the instructions specified in the manual.

FAQs

Q: Which one is better in terms of weight reduction, liposuction or this equipment?

A: The principle of liposuction is sucking out excessive fat from a certain body part through vacuum suction so as to achieve the purpose of making local parts lose weight quickly. It has an instant effect, but requires anesthesia during the operation, and has recovery period. It's an operation and there is risk. However, MS-12C2 is zero-risk and without side effects, and can strongly blast fat and effectively reduce local fat. Women can become thin while at home. And it can make women more sexy and enhance female charm.

Q: Will I experience rebound after the treatment?

A: MS-12C2 is an equipment that can strongly blast fat. It will not rebound that easy after reduction. Since ultrasonic fat-dissolving is removing fat instead of water, and it takes a long time for the fat to be accumulated, thus the weight will not get rebound easily.

Q: Do I need to be on a diet?

A: You need to keep on a diet to some extent since ultrasonic treatment is followed by fat blasting and accelerated metabolism. Avoid eating spicy, oily, or fried foods so as not to affect metabolism. Proper exercise can be done for sweating, and make it discharge successfully from your body. In this way the effect of weight loss will be more obvious.

Q: Why does drumming in the ears occur?

A: It's because ultrasonic wave has a strong vibration and with a frequency above 60KHZ. It acts on deep subcutaneous fat layer through focusing. Ultrasonic wave with concentrated energy effectively causes high speed friction among fat cells in focal area, resulting in bursting, heating and eventually in breaking and emulsifying, during which there is sound. Therefore, slight tinnitus will accompany during the treatment.

Operational Diagrams


Waist & Abdomen Fat-dissolving & Sculpting




1. Wet the towel with warm water, then clean the parts to be treated, 2 minutes.




2. Apply gel evenly to the parts to be treated, 1 minute.

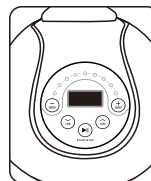


3. Set energy level between 2 and 7 and time between 15 to 20 minutes, and adjust them according to individual tolerance, 1 minute.

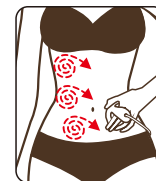


4. Move in small circles and evenly on parts with fat, 3 minutes.

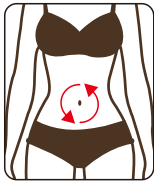

5. Move in big circles around belly button, 3 minutes.

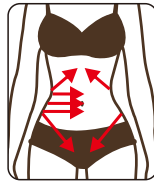

6. Push from the side of waist to the middle of abdomen, 3 minutes.

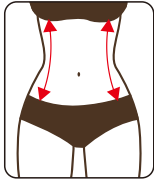


7. Slide up and down on the two sides of waist, 3 minutes.



8. Clean the treated parts, 1 minute. Treatment done.

Arms Fat-dissolving & Sculpting

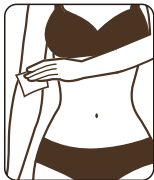


1. Wet the towel with warm water, then clean the parts to be treated, 2 minutes.




2. Apply gel evenly to the parts to be treated, 1 minute.




3. Set energy level between 2 and 7 and time between 15 to 20 minutes, and adjust them according to individual tolerance, 1 minute.

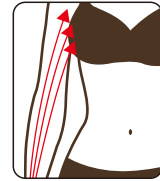

4. Push from lower arm with too much fat to armpit along the three inside and outside channels and collaterals, 5 minutes.

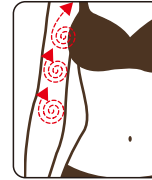

5. Move in small circles on lower arm with fat, 3 minutes.

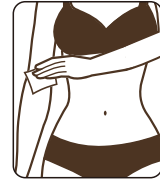

7. Clean the treated parts, 1 minute. Treatment done.

6. Do the same with the other arm.

Legs Fat-dissolving & Sculpting




1. Wet the towel with warm water, then clean the parts to be treated, 2 minutes.




2. Apply gel evenly to the parts to be treated, 1 minute.




3. Set energy level between 2 and 7 and time between 15 to 20 minutes, and adjust them according to individual tolerance, 1 minute.

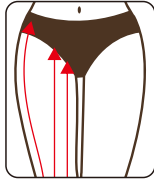


4. Push from leg part with fat to thigh root along the three channels and collaterals, 3 minutes.

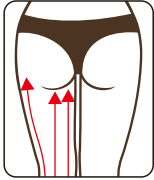

5. Push from leg part with fat to thigh root along the three channels and collaterals. 3 minutes.

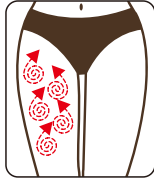

6. Move in small circles on parts with fat, 3 minutes.

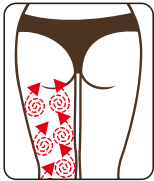

7. Move in small circles on parts with fat, 3 minutes.

8. Do the same with the other leg.



9. Clean the treated parts, 1 minute. Treatment done.

Course of Treatment Recommended:

For initial stage, the treatment can be done for 3 to 4 times in a week. Later, it can be done for 1 to 2 times in a week to consolidate the effect. One course consists of ten treatments. After one treatment, the fat reduces obviously, and flab of the treated part lessens to some extent. After two treatments, the fat reduces even more. Treatment can be chosen as per individual demands.

This equipment can be used to treat arms, abdomen and waist, and legs.