


# SMART V

JUST FOR YOU

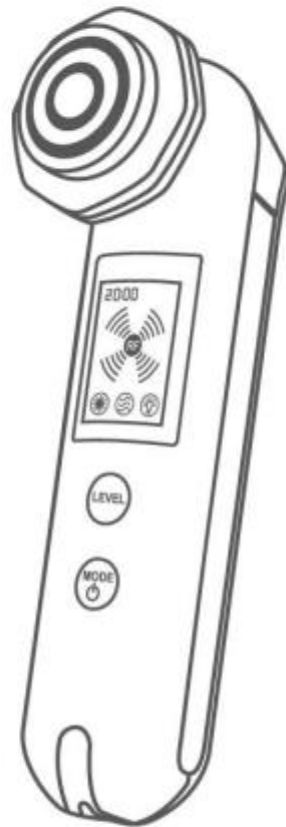

## **SMART V Skin Renewal System** **RF/LED Photon/Massage** **Instruction Manual**

Thank you very much for choosing this product-SMART V

Please read this manual carefully before use

This product is intended for home use and cosmetic purpose only, using this product for treating any medical condition or illness is not supported by the manufacture.

# SMART V Skin Renewal System

Congratulations on your purchase

To fully benefit from the support that our offers, please read the user manual carefully before you use the appliance and save it for future reference.

Thank you for purchasing our new SMART V device-3IN1 Skin Renewal System. This new clinical skin renewal device is adapter for home use, so you can enjoy professional, clinical, aesthetic skin treatment results at home with a safe, effective and easy-to-use device. For the first time, you'll experience real gain without any pain.

Our SMART V is the only scientifically proven home device that generates a real change, starting from within. Treatments will tighten the skin, reduce fine lines and wrinkles, leaving it radiant, youthful and glowing.

Powered by the revolutionary patented SMART RF technology that focuses RF energy deep into the dermis, much like a laser lens focused lights. The SMART RF power stimulates dermal activity, tightening collagen fibers and increasing new collagen production. The skin becomes tighter, more supple, with significant reduction of fine lines and wrinkles. In a short time, it effectively reverses the signs of aging, easily adjusting to all skin types and colors.

- **The best alternative to face-lift surgery, botox and other professional procedures.**
- **Painless, safe and pleasant, offering an enjoyable self-treatment at home.**
- **Scientifically and clinically proven to produce real, measurable and long-lasting results.**

1. Safety Caution .....	02
2. Core Technology .....	02
3. Parts .....	03
4. Function Introduction .....	03
What It Does.....	03
Principle of Instrument.....	04
• SMART RF Therapy.....	04
• Photon Therapy .....	04
• Vibration Massage .....	04
5. How TO Use .....	04
What to Prepare.....	04
Introduction to the keys .....	05
Daily Care .....	05
Skills .....	08
Expected Results .....	11
Care and Maintenance .....	11
Specification.....	11

# 1. Safety Caution

This new clinical SMART V skin renewal device is adapted for home use. Please use it for intended purpose, as described in this user manual only.

## Warnings

- Never leave device unattended when plugged in.
- Do not use the device if it is disassembled.
- Use the device only with the supplied cord and power supply.
- Keep out of reach of children.
- Important! Do not use the device over eyelids, lips, ears, or genital areas.
- Don't use it if any of the following apply to you.
- Under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Current state or history of cancer, especially skin cancer or premalignant moles.
- Any active skin condition in the treatment area, such as sores, psoriasis, eczema, or rash.
- If you have side effects in the treatment area, please consult with a doctor as soon as a possible.


**Note:** In any case of discomfort or uncertain potential side effects, please consult with a doctor

## 2. Core Technology


- Focuses SMART RF power deep into the dermis to generate heat and transition from within.
- Uniquely designed SMART RF electronic tip focuses power deep into the dermis, much like a less focused light.
- This increased dermal activity results in visible skin tightening, firming and rejuvenation. firm lines and wrinkles are significantly reduced and the skin becomes smoother.
- **Three times per week, for a period of 6-8 weeks to begin seeing visible skin improvements.**
- **To Achieve long-lasting results, continue repeating the treatment twice per week.**
- We commit to safety, ease of use and high efficiency principle. It is important to read and follow all instructions and safety precautions described in this user manual .
- Designed for all skin types and colors.

**Note: Please read this manual carefully before use!**


### 3.Parts

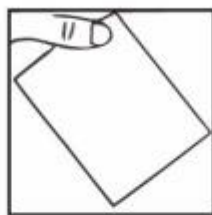

cotton pad

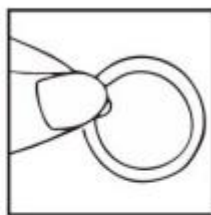


cotton holder

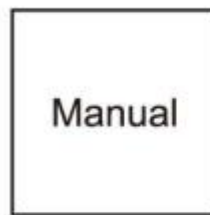

Manual

Manual

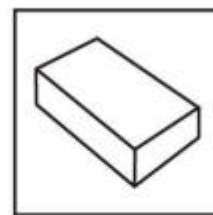

Gift Box

### 4.Function Introduction

#### What is does

This instrument focused on the following problems.

- Tighten sagging skin, brighten skin tone.
- Reduce wrinkles and fine lines, smooth cellulite.
- Improve skin elasticity and skin texture.

## Principle of Instrument

### **SMART RF Therapy**

It utilizes a special technology of touch sensor, can commend inner RF procedure to emit RF energy from low to high adjusted automatically upon skin's temperature, focuses power which can up to 40-50 celsius deep into the dermis, thus dermis collagen may immediately shrink after wrinkle treatment, client may immediately sense the skin tightening effects as it is being lifted and firmed.

### **Photon Therapy**

Red LED photon with wavelength of 625 $\pm$ 3nm can whiten and lighten dark spots, tender skin and dispel wrinkles, repair damaged skin, heal slender wrinkles, tighten pores, and proliferate collagen protein. Blue LED photon with wavelength of 465 $\pm$ 3nm can diminish inflammation, battle bacterium and propionibacterium acne, treats acne and scars.

### **Vibration Massage**

Vibration massage can relieve acne, relieve fatigue, and refine skin molecules with 3 million hertz to help deliver skin care nutrients deep into the skin, effectively tighten skin.

## 5. How to use

### **What to prepare**

#### **Personal Preparation**

The best time for treatment is after bath and before sleep. Generally, the ideal results can be achieved after 2 weeks of treatment following the recommended frequency. Later, you can reduce the frequency as according to your skin condition.

**Test this equipment on the wrist before using on your face.**

#### **Device preparation**



Firstly, Connect the device with power supply correctly, make sure it is well charged.





Secondly, please take one piece of cotton on the probe, then fix the cotton by cotton holder.



Thirdly, Press the ON/OFF button for about 4 seconds to switch the device.

## Introduction to the keys



**Power/Mode Power on/off**, it is connected with the power supply to be charged. Mode selection: Press "MODE" to select the desired function, the function will be flashing when selected; It will automatically shut down in 6 seconds if there is no further operation.



**LEVEL** Level adjustment: Select "MODE", Press "LEVEL" key to adjust the levels RF- 5 Levels.



88:88 — Working time (default 30 minutes)

RF — RF: 5 Levels

— Vibration Massage

— Photon: Red Blue two modes

## Daily Skincare

- Here is just a general introduction of home use, for more details, please consult a professional beauty consultant or a doctor.
- Apply 4-5g cosmetic essence such as model OT-AC073 on each treatment area.



**Tips:** Please keep the SMART V probe in full contact with the skin during treatment, move the probe on the surface of the cared area, in an even and thorough way.

**Important:** The most important thing is to choose suitable cosmetics for your skin. It is not only for therapy, but also an important foundation to provide nutrient to your skin. Therefore, be careful to select the right cosmetics.

## Daily Face Care

Use three times per week at most, no more than 30 minutes each time to achieve the best result. Choose function & levels according to your skin type.

### **Step1: Skin clean**

Remove make up, wash face, facial tone and then facial spray to open skin pores.

### **Step2: Nutrition input(create a perfect feature)**

Turn on the machine, apply 4-5 g cosmetic essence on the cotton pad, fix the cotton pad by cotton holder, select the "RF" function, then select a comfortable level. It is better to use water liquid cosmetic (toner, moisturizer, water, essence, essential oil, aloe vera gel or other liquid). For best nutrition absorption, move slowly from the chin to your ears for 10-15 seconds, then move from left to right on the forehead area.

### **Step3: Photon rejuvenation therapy**

Select "Photon" and choose red light, use a point-by-point moving method, especially for dark spots and wrinkles; For red light irradiation move every 10-15 seconds. Photon rejuvenation care session is recommended to be 10-15 minutes.

### **Step4: Vibration massage**

Choose "Massage" function, select "LEVEL", use the point-by-point moving method for about 10-15 minutes.

### **Step5: Basic skin care**

Refer to the picture as following "Skills".

## How to work best Tighten sagging skin(every3 days)

### **Step1: Deep clean the skin**

Remove make up, wash face, facial tone and then do facial spray to open skin pores.

**Step2: Nutrition for anti-aging importing, skin beautification ( same as you do your daily care)RF Level:4-5**

### **Step3: Import nutrition gel/mask**

Select Red LED photon use a point-by point moving method, especially for sagging treatment area, work time should last for about 10-15 minutes.

**Step4: Moisturizing mask working time:15 minutes**

### **Step5: Vibration massage**

Choose "Massage" function, select "LEVEL", select Blue LED photon, use the point-by-point moving method for about 10-15 minutes.

### **Step6: Basic skin care**

Refer to the picture as following " Skills".

## Acne Care (Once every night )

### **Step1: Deep clean the skin**

Remove make up, wash face, facial tone and then do facial spray to open skin pores.

**Step2: Nutrition for acne importing, skin beautification ( same as you do your daily care) RF Level: 2-3**

### **Step3: Import nutrition gel/mask**

Select Blue LED photon use a point-by point moving method, especially for sagging treatment area, work time should last for about 10-15 minutes.

**Step4: Acne care mask working time:15 minutes**

### **Step5: Vibration massage**

Choose "Massage" function, select "LEVEL", choose Red LED photon , use the point-by-point moving method for about 10-15 minutes.

### **Step6: Basic skin care**

Refer to the picture as following " Skills".



# Eye edema care(weekly)

## Step1: Deep clean the skin

Remove make up, wash face, facial tone and then do facial spray to open skin pores.

## Step2: Nutrition for eye edema importing, skin beautification ( same as you do your daily care) RF Level:1-2

## Step3: Import nutrition gel/mask for eye edema

Close your eyes, select Red LED photon use a point-by point moving method, especially for eye edema treatment area, work time should last for about 10 minutes.

## Step4: Tighten mask for eyes working time: 15 minutes

## Step5: Vibration massage

Choose "Massage" function, select "LEVEL", choose Blue LED photon use the point-by-point moving method for about 10 minutes.

## Step6: Basic skin care

Refer to the picture as following " Skills".

## Skills



### 1.Neck

Move in small circular motions from the inside of the neck to the outside, then lift from the neck upward to the cheek, neck treatment area can be divided into three lines(up,medium,down areas) to operate.

Attention: strength should be top heavier and down lighter, operation electrodes should tight skin, repeat 8-10 times.



## 2.Chin

Make operation head move in small circular motions by pressing tightly against the double chin and then pull back to the back of the ear, repeat 5-8 times.



Make operation head move in small circular motions by clinging to the chin and circle out to the ear beads. Then pulled down to the chin with light force,repeat 5-8 times.



## 3.Mouth

Make operation head move in small circular motion by clinging to the corner of mouth, and circle out to the ear beads. Then pulled down to the chin with light force,repeat 5-8 times.



## 4.Cheek

Make operation head move in small circular motion by closing to the light bone of nose, circle out to the ear beads. Then pulled down to the chin with light force, repeat 5-8 times.



### 5.Forehead

Make operation head move in small circular motion by closing to the light bone of nose, circle out to the ear beads. Then pulled down to the chin with light force, repeat 5-8 times.



### 6.Eyes' around

Dredge the meridian and collaterals of eyes, repeat 5-8 times.



Move in small circular motion in the lower eye then lift to the temple, repeat 5-8 times.



Move in small circular motion from the corner of the eye to the Jingming points to the eyebrow, repeat 5-8 times.

## Expected Results

A small amount of erythema and slight redness may occur on the skin area after treatment. This is quite normal and phenomenon will disappear shortly. Skin will be obviously more elastic after a session of treatments combined with proper cosmetic. The wrinkles and fine lines will be improved.

Clinical research and scientific measurements show that after 6-8 weeks of use, 30% of fine lines 16% of wrinkles can be reduced. Specific effects vary based on individual skin characteristics. If you follow the regiment and continue doing it, the skin will show amazing results.

## Care and maintenance

- Before any cleaning attempts, make sure the power is off.
- Disconnect the device from the power supply before cleaning it.
- This care instrument is suitable for private use only, do not share it with another individual.
- Keep the system in a cool, dry place.
- Clean the residue of gels with a soft damp cloth or paper towel immediately after each use,
- **Alcohol is prohibited to clean or wipe the device.**
- Keep the device dry & clean at all times.

## If the product does not work properly, check the following

- Ensure the power cord is well connected with the device.
- Confirm LCD screen light up.
- Select function according to the instructions in this manual.
- If the power cord is damaged, contact your vendor for a replacement.
- If you did not select any operation within 6 seconds, the unit will automatically shut down to protect the circuitry.
- Ensure the working time, the default working time is 20 minutes.

## Specification

Model: SMART V

Net Weight:

Size: 6.6" \* 1.7" \* 1.2"

Working Environment: 5-40°C

Power supply: 100-2240V 50/60Hz

Input: AC100-240V 50/60HZ 0.6A

Output: DC9V 2A

US, AU, UK, EU adapters offered

Warranty: 1Year