

User Manual

SR-9000



Preface

Dear users:

We're pleased to present to you the *Galvanic Current Therapy Personal Use Micro-current BIO Beauty Apparatus* that mainly focus on dry skin, clogged pore, enlarged pore, inelastic skin, etc. It aims to deliver a safe and effective beauty experience that helps to replenish skin with nutrition, extract skin gunks and debris through high-tech means. *Galvanic Current Therapy Personal Use Micro-current BIO Beauty Apparatus* is specialized skin care machine for professional use, which needs to be handled by professionals with proper training. Any improper use will result in adverse consequences. Therefore we advise any personnel to read this guide thoroughly and follow the instructions strictly before operation.

We believe our quality product will bring in good returns on your end. Also, you can count on our worry-free after-sales service.

Thank you!

Table of Contents

Preface
Part I
Brief Introduction
Advantages
Working Principles
Main Effects
Applicable Range
Inapplicable Range
Notes After Treatment
Part II
1. Detailed Operations
2. Technical Parameters
3. Safety Precautions
4. Dos and Don'ts
5. Troubleshooting & Solutions
6. FAQs
7. Packing List
8. Operational Diagrams

Part I

Brief Introduction

Galvanic Current Therapy Personal Use Micro-current BIO Beauty Apparatus is currently the most popular beauty gadgets that work on skin cleaning, hydrating, tightening, and anti-wrinkle. Ergonomic design brings you a comfortable, effort-saving experience. It bring out skin debris and lead in nutrition. With its powerful effect, it can replace traditional cleansing, facial skin care programs at SPA. Moreover, it's convenient and easy to operate. It solves various skin issues for beauty-lovers. This device requires no injection, no medication, no operation, thus having no side-effects. It is used externally throughout the entire process and comes with immediate result. Ions use direct current polarity to strengthen the penetration of nutritive active substances into the skin, so as to beautify the skin, and separate the dust, garbage and excess cutin in the skin for deep cleaning of the skin.

Advantages

1. Small, convenient, portable and easy to operate.
2. Painless, non-invasive during the process. Skipping any recovery periods without affecting any work and life plans.
3. Simple operation, easy to learn, and remarkable effect.
4. Two operational heads to choose from based on skin needs.
5. No unevenness, bleeding, swelling or stasis.
6. It uses the interaction between positive and negative ions to produce mutual attraction and mutual repulsion, achieving beauty results of deep cleaning and leading in nutrition.

Working Principles

Ion+/- nutrition lead-in uses the principle that like charges repel each other; unlike charges attract so that nutrient active substances of the same polarity can penetrate deep into the skin through the sweat gland or help strengthen the penetration so as to achieve the effect of purification of the skin.

ION- : When it's at ION- mode, it can be used for nutrition lead-in. The customer holds the '+' passive electrode probe while the probe is at it's '-' mode so the nutrients applied on the skin is negatively charged. When the active substances in essential oil/cream or other skin care products meet the charges of same polarity will attract and therefore accelerate the process of penetrating through pores, sweat gland, and eventually fasten the nutrients absorption. In this way, it can achieve it's skin care effect.

ION+ : When it's at ION+ mode, it can be used for waste lead-out. The customer holds the '-' passive electrode probe so some of the dust, dirt and excess cutin in our skin is negatively charged. Through the equivalent exchange of ions (positive for positive, negative for negative), the gunks hidden in our pores can be easily led out, thus achieving the effect of cleaning. In this way, it can achieve its skin cleansing effect.

Main Effects

1. Deep cleanse skin and exfoliate.
2. Boost skin absorption and deep hydrate.
3. Even out fine lines and fade wrinkles.
4. Improve uneven skin tone and dull skin.
5. Improve soft and inelastic skin.
6. Lead out gunks hidden inside pores and minimize pores.
7. Improve coarse skin and refine skin.

Applicable Range

1. Those with dry, dehydrated and malnourished skin.
2. Those with poor skin absorption and lack of moisture at skin base.
3. Those with dull skin and uneven skin tone.
4. Those with clogged pores and thick stratum corneum layer.
5. Those with soft and inelastic skin.
6. Those with coarse skin and enlarged pores.

Inapplicable Range

1. Those with "3 High" (hypertension, hyperlipidemia, hyperglycemia) and heart disease.
2. Those who have pace makers and metal implants.
3. Those who are in their pregnancy, menstruation or lactation.
4. Those who have skin disease and infectious disease.
5. Those who are allergic to electric current.
6. Those who have extensive facial trauma.
7. Those who have just received injection, such as micro-needling, intradermal injection, micro plastic surgery or those who are still recovery from surgeries.

Notes after Treatment

1. Strengthen hydration and always use sun protection.
2. Do not use stimulating skin care products including scrub, exfoliant and AHAs within one week after operation.
3. Avoid staying up late, smoking, drinking, and eat more fruits and vegetables and less greasy food.
4. Apply a hydrating sheet mask at least 3 times a week.
5. Do not wash your face with hot water within 3 days after operation. (Use warm water instead)

Part II


1. Detailed Operation

Detailed installation is as follows:





Detailed Functions

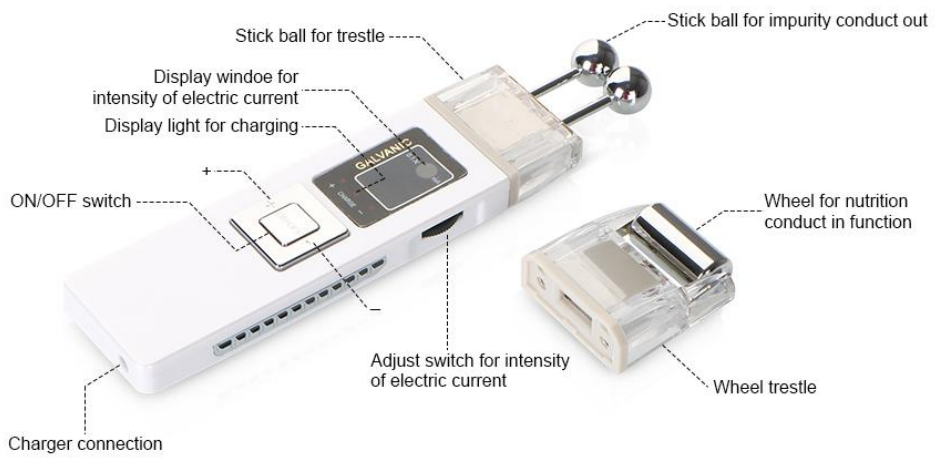




ION + Probe

It can be used to clean skin, exfoliate, bring out debris, dirt and can be used together with exfoliating products.



ION - Probe

It can be used to help skin absorb nutrients, tighten skin and shrink pores. It can be used together with skin care products, e.g. essence, lotion, nutrients, etc. Pair with different product based on skin needs.

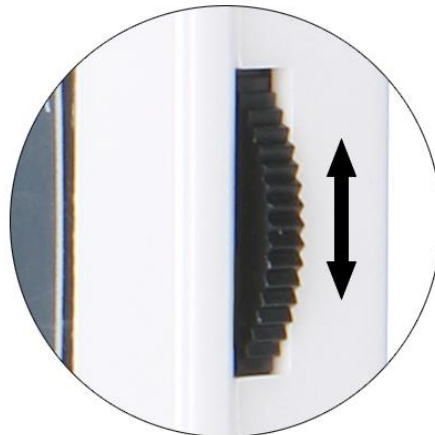

ION Intensity +/-

Slide up to decrease ION intensity and till the end to turn the ION level into none.
Slide down to increase ION intensity and till the end to turn the ION level to the

maximum.

ION Intensity Range: 0.5-10

Reference Level: 2-8

Higher level comes stronger electric feeling. Adjust level based on one's own comfort level.



ON/OFF

Click the ON/OFF button to turn on/off the device. Select “-” to lead out gunks and lead in nutrition. “+” mode brings less felt electric feeling and less pronounced results.



Make sure your hand is in contact with the metal piece on two sides and the probe in in full contact with the skin



The device cannot work when your hand makes no contact with the metal on two sides



The device doesn't work until the probe is in full contact with skin



Normal Working State

Red or blue light turns on after the probe makes full contact with the skin

2. Technical Parameters

Voltage: 100-240 V (Worldwide useable)

110V(USA & Canada) & 220V (Europe & U.K.&AU)

With chargeable lithium battery, can use easily at home or during travel.

Frequency: 50-60Hz

Galvanic Current level (Power level): 0.05mA--1.0 mA (Adjustable)

Color: Silver White

3. Safety Precautions

Please consult your doctor or professionals before using the device if you have the following symptoms:

1. Those who are in their pregnancy, menstruation, or lactation.
2. Those who have pace makers, metal or other electrical appliances implants.
3. Those who have epilepsy, heart disease, severe "3 High" (hypertension, hyperlipidemia, hyperglycemia), diabetes, hyperthyroidism and other severe diseases.
4. Those who have skin diseases and infectious diseases.

5. Use with caution for those who are allergic to electric current.

4. Dos and Dont's

1. Ensure that the voltage of the device is appropriate. If the voltage of the local power supply is unstable, we recommend that users increase to the matching power voltage.
2. In order to ensure the therapeutic effect and normal service life expectancy of the device, please use the specified accessories provided or recommended by the original manufacturer.
3. The device should not be placed in wet environment or near water, nor should they be directly exposed to sunlight.
4. Do not place the device near a strong heat source, because this may affect the life and normal use of the device.
5. Please restrain from using the device on eyes, thyroid, parathyroid, testicles, abdomen of pregnant women, and pacemakers.
6. Patients who are suffering from illness should consult with a doctor and gain permission from the doctor before using.
7. Use device or train device operators in strict accordance with instructions in the user manual.
8. Start from the lowest energy level and slowly add up.
9. When using this device, the operating parts must be kept moist and dry skin should be avoided.
10. Charge in time when the battery is low. Avoid half treatment to effect the curative result.
11. Use with caution for those who are allergic to electric current.
12. Do not soak the device into any kind of liquid.
13. Make sure hand is in contact with the metal on two sides of the device.

5. Troubleshooting & Solutions

1. I can't turn it on?

- A. Please check if the device is fully charged.
- B. If the above methods cannot address the issue, please contact the instrument dealer for assistance.

2. I can't feel much during the treatment after turning on the device?

- A. Please check whether the probe comes loose or not installed properly.
- B. Please check if the energy level is set at the lowest level and readjust.
- C. Please make sure there's media, e.g essence, toner, lotion, etc.
- D. Please make sure the probe is in full contact with skin.
- E. Please make sure hand is in contact with the metal on two sides of the device.

F. Please make sure you select “-” mode.

G. If the above methods cannot address the issue, please contact the instrument dealer for assistance.

6.FAQs

1. Q: Is ion treatment harmful to the skin?

A: No. Because ion+/- nutrition lead-in uses the principle that like charges repel each other; unlike charges attract so that nutrient active substances of the same polarity can penetrate deep into the skin through the sweat gland or help strengthen the penetration so as to achieve the effect of purification of the skin. So ion treatment does not harm the skin. On the contrary, it's a great approach to skin care.

2. Q: Can I wear makeup after the treatment.

A: Yes. The electric ion causes the chaos cell at the inner layer of facial skin to rearrange through the positive ion so it precipitates the toxin, removes the skin inner layer dirt, so the skin can breathe again. It's non-invasive, no wound, safe, healthy, fast and effective.

3. Q: Will my skin deteriorate if I stop doing it?

A: No. This only targets at certain skin problems so that it can improve and deeply clean the skin, which belongs to the maintenance and improving skin of malnutrition caused by other factors. It helps to absorb nutrients in skin cares products and relieve sallow, dehydrated, coarse or wrinkled skin. Your skin will only go back to its natural state if your stopped the treatment.

7. Packing List





1 x Main Machine

1 x Ion + Probe





1 x Ion - Probe

1x Power Cord

8. Operational Diagrams

Parameter Adjustment	Product	Techniques	Diagrams
Cleaning&Hydrating: 60 min, once/week			
<p>Probe: ION- probe ION+ probe</p> <p>Mode: “_”</p> <p>Energy Level: 2-8</p>	<p>Makeup remover+ Cleanser+ Exfoliating cream+ Hydrating essence+ Sheet mask</p>	<ol style="list-style-type: none"> 1. Remove makeup thoroughly. 2. Clean face with cleanser. 3. Apply exfoliating cream all over the face. 4. Select ION+ probe and install. 5. Adjust mode and energy level. 6. When you hold the device, make sure your hand makes contact with the metal on two sides of the device. Start the treatment. 7. ION+ Cleaning Move the probe slowly from chin to earlobe, mouth corner to ear gate, nose wing to temple. Slide back and forth, 3-5 times. 8. Move slowly back and forth along the nose wing, 3-5 times. 9. Do the same on the other side. 10. Slide back and forth from forehead to hairline, 3-5 times. 11. Pause the device and switch to ION- probe. 12. Wash face clean with lukewarm water. 13. Apply hydrating essence 	<p>Technique 7</p>  <hr/> <p>Technique 10</p>  <hr/> <p>Technique 16</p>  <hr/> <p>Technique 19</p> 

		<p>evenly all over the face, 1 min.</p> <p>14. Readjust mode and energy level.</p> <p>15. When you hold the device, make sure your hand makes contact with the metal on two sides of the device. Start the treatment.</p> <p>16. ION- Hydration</p> <p>Move the probe from chin to earlobe, mouth corner to ear gate, nose wing to temple in parallel lines to lift the treatment site, 3-5 times.</p> <p>17. Slide slowly back and forth from nose wing to inner eye corner, 3-5 times.</p> <p>18. Do the same on the other side.</p> <p>19. Move toward hairline on forehead to lift the treatment site, 3 times.</p> <p>20. Slide slowly back and forth from left to right, 3 times.</p> <p>21. Apply sheet mask, 15 min.</p> <p>22. Remove the mask and wash face clean with lukewarm water.</p> <p>23. Apply toner, lotion, cream and sunscreen.</p>	
<p>Recommended treatment:</p> <p>It is advised to start with 3-4 times a week. After 2-3 weeks, switch to 2-3 times per week. One session helps clear pores and hydrate skin. One month of persistent use helps enhance the effect. After 2 months, skin appears to be hydrated, brightened and elastic. Three months of treatment helps refine skin. Long term use helps skin maintain glossy and clear and prevent congestion.</p>			
<p>Skin Firming & Anti-wrinkle: 60 min, once/week</p>			
<p>Probe: ION+ probe</p>	<p>Makeup remover+</p>	<p>1. Remove makeup thoroughly.</p>	<p>Technique 7</p>

<p>Mode: “_”</p> <p>Energy Level: 2-8</p>	<p>Cleanser+ Sheet mask+ Anti-wrinkle essence</p>	<ol style="list-style-type: none"> 2. Clean face with cleanser. 3. Apply anti-wrinkle essence evenly all over the face. 4. Select ION- probe and install. 5. Adjust mode and energy level. 6. When you hold the device, make sure your hand makes contact with the metal on two sides of the device. Start the treatment. 7. ION- Anti-Wrinkle Slide back and forth along the jawline, 3-5 times. 8. Move the probe from chin to earlobe, mouth corner to ear gate, nose wing to temple in parallel lines to lift the treatment site, 3-5 times. 9. Slide from the inner eye lid to temple, tracing the lower eyelid, 3 times. 10. Lift from the eye corner towards hairline, 3 times. 11. Do the same on the other side. 12. Move towards hairline on forehead to lift the treatment site, 3-5 times. 13. Apply sheet mask, 15 min. 14. Remove the mask and wash face clean with lukewarm water. 15. Apply toner, essence, eye cream, cream and sunscreen. 	<div style="text-align: center;">  </div> <p style="text-align: center;">Technique 8</p> <div style="text-align: center;">  </div> <p style="text-align: center;">Technique 9</p> <div style="text-align: center;">  </div> <p style="text-align: center;">Technique 12</p> <div style="text-align: center;">  </div>
---	---	--	--

Recommended treatment:

It is advised to start with 3-4 times a week. After 2-3 weeks, switch to 2-3 times per week. No visible results may be seen after one session. A month of persistent use helps brighten skin, shrink pores and skin appears to be translucent. Two months of treatment help refine skin, decrease fine lines around eyes, fade wrinkles and increase skin elasticity. Three months of use help even out skin tone, and tighten skin. Long term use helps prevent loose, sagging skin and delay skin aging.

