



WL-7005C

# Instructions

## Warning

Dear users:

Thanks for choosing our company's product, in order to fully make use of the products, we sincerely suggest follow:

- ◆◆ Please read the instruction in detail and keep it well for further reading and consulting.
- ◆◆ Please follow the instruction to fit and operate the machine correctly.
- ◆◆ Please do not remove and change any accessories of the machine.
- ◆◆ It forbids others to open and remove the mainframe except the authorized people in the company.
- ◆◆ It forbids others to change and maintain the machine except the professional people.
- ◆◆ Please contact us if there is any trouble during operation.

### Please read the NOTES carefully before use Machine.

1. Please make sure the tank is filled with water before start machine .
2. Please put at least two pieces membranes on treat area before start treatment .
3. Please clear the filter in time .otherwise it will cause residual liquid to be poured into the machine and burned out parts .
4. we recommend each area treat time do not over 30mins.over 30mins ,please do not allow clients stay alone,need to check client situation in time.
5. During treatment ,check skin condition every five minutes.
6. If can not feel cold in one minute , please stop and check the condition of machine.
5. Please change water once a day.

6. Patients using the following medical electrical equipment can not use the equipment.

- Implanted cardiac pacemakers or other implanted electronic devices;
- Life-sustaining artificial heart-lung machines and other equipment;
- Portable ECG measurement equipment;

The following patients should be agreed by the professional doctors before using the equipment.

- Patients with acute
- Malignancy
- Patients with Infectious
- Menstrual period women
- Pregnant women
- Patients with heart
- Patients with fever
- Patients are receiving operation or those with abnormal body phenomenon;

After Sale service

◆ We offer ONE YEAR warranty. Within warranty period, all exchanges or replacement due to defects are free. (Buyer is responsible for shipping charges).

If the machine is over warranty, please also do not hesitate to contact us. we will assure you of our best service all the time.

◆ How to contact us?

Phone:(+86) 189 -487-55929

Tel: 0755-86368500

Email: sales@mychway.com

This product does not belong to medical, please don't use for disease treatment!

## **Chapter 1: Safety Warnings**

### **I. Electricity & Machinery Safety**

1. The machine is universal general voltage. Can work with 100V & 200V voltage.
2. Please keep the voltage at stable condition.

### **II. Contraindication Safety**

1. People with heart disease, high blood pressure, or configured cardiac Pacemaker.
2. People with acute inflammation, asthma, deep vein thrombosis, thyroncus, cancer.
3. People with hemorrhagic diseases, trauma or who is bleeding.
4. Pregnant women
5. Medical plastic artificial or metal parts inside the body such as metal teeth or silicone breast implants.
6. People with abnormal immune system.
7. People with numb or insensitive to heat.
8. If you are feeling unwell, Prohibited in the ears, nose, eye throat and do not use this device in advance and tell your doctor.
9. Children under 12 years old are not allowed.

### **III. Treatment Safety**

1. Treatment Time is less than 30mins/Body Part/Each Treatment.
2. The energy setting is from low to high gradually which should be subject to skin comfortableness for different people.
3. The patient had better wear the protective eye glasses.
4. If you cannot tolerant heat or uncomfortable, please stop the treatment in advance and turn to your doctor for help.

## **Chapter 2: General Introduction**

### **Fat Dissolve Vacuum**

Why the Cool Fat Dissolve Vacuum Procedure Is Different?

It is the non-invasive cooling of adipose tissue to induce Fat Dissolve - the breaking down of fat cells - to reduce body fat without damage to other tissues.

It is fundamentally different from other non- or minimally invasive modalities. Other methods of

fat removal primarily involve necrotic cell death by damaging fat with heat, high-intensity focused ultrasound, or chemical injections. Each approach poses potential technical challenges, particularly with respect to targeting the right tissue depth and unintended damage to other structures close to or within the fat layer. Predictability of efficacy using these other techniques is also uncertain. In contrast, the Cool Fat Dissolve operation induces apoptosis only in fat cells to gently and gradually reduce the fat layer while preserving all other tissue.

As triglyceride in fats will be converted into solid in particular low temperatures, it uses advanced cooling technology to selectively target fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues, reduce unwanted fat. When fat cells are exposed to precise cooling, they trigger a process of natural removal that gradually reduces the thickness of the fat layer. And the fat cells in the treated area are gently eliminated through the body's normal metabolism process, to eliminate unwanted fat.

#### **40K cavitation**

What is the Ultrasound Cavitation operation?

It uses a 40kHz focused ultrasound cavitation effect, to cause tons of microscopic air bubbles to surround the LEDcyte (fat cell) membranes and with its resonant frequency, the bubbles cause a high enough to cause the fat cells to break up instantly.

The broken LEDcytes are then moved to the liver and safely metabolized out of human body which is assisted by drinking a lot of water.

It produces fatty acids from adipocytes (fat-laden cells) due to a phenomenon known as stable cavitation, and increased cell permeability of the adipocytes. You get an increased flow of triglycerides. Furthermore, there is a breakdown of the fibers through mechanical action (jet stream). Finally, some content of adipocytes is metabolized by the liver and eliminated in the urine.

#### **Benefits Of RF Facial/Body Treatment**

Pure RF energy is used to treat skin laxity by facial tightening, and is best suited to patients with mild to moderate sagging of facial tissues, usually those in their mid thirties to fifties.

It is commonly used to treat the forehead area (lift brows and tightens forehead skin), under the eyes (thickens skin and tightens and lifts the upper cheek skin), cheeks (lifts slightly and diminishes enlarged pores), mid-face, jaw line (diminishes jowling and sagging along the jaw line), and the neck (tightens skin and diminishes wrinkling of the neck).

#### **You can get 5 benefits from this item:**

1. Fat Reduction
2. Body remodeling
3. Cellulite operation and skin partial reconstruction
4. Increase Collagen Level
5. Tightening loose skin

### **How does 650NM LIGHT system work?**

It involves the application of a highly sensitive led to dissolve the fat deposits situated in the upper layers of the human skin. The energy emanating from the led breaks up the fat cells to produce an oily material.

It penetrate into skin surface and stimulates the fat cell membranes,changing their permeability,the fat cells reduce their overall size and intracellular fat is released then the fatty trigly cerides flow out of the disrupted cell membranes and into the interstitial space,where they gradually passthrough the body's natural metabolic functions with no harmful physiological effects,this process is resulting in inch loss of patients.

### **Technical Specification**

Rated input voltage:AC220V-240V 50HZ /AC110V-120V 60Hz

Rated input power:350VA

Cool cold slimming Vacuum

Power:Up to 200W

Cooling device output temperature:5~ -5° C

Cooling device output :0-5Kpa

Cooling liquid: pure water

Cavitation frequency : 40KHz

Cavitation Power:up to 60W

Cavitation tip:60mm diameter

RF for BODY

RF frequency:3MHz

RF power:up to 30W

RF for EYE

RF frequency:3MHz

RF power:up to 30W

Lipo laser wavelength :635nm

Each diode:160mw

laser Pads: 12 with 8 Large(8 Diodes/Pad)+4 Small(2 Diodes/Pad)

## Packing Included

- 1 x Host machine
- 1 x Cool Vacuum Head
- 1 x 40KHz Cavitation Head
- 1 x Sixpolar RF head for body
- 8 x Big Head(each with 10 light)
- 4 x Small Head(each with 2 light)
- 3 x Holders for All parts

## Chapter 3: Preparation before use


1. Install the attached on-hook on the machine, put the probes into the corresponding racks.
2. Before using the Cryolipolysis Cool-sculpting Slimming Machine, it is required to connect the funnel with inlet hole and open the cap of overflow hole before adding distilled water or purified water.
3. Fill machine with water by funnel through the inlet until water out from the overflow continuously, and then take off the funnel and close cap with the overflow hole.




4. Check the host power cable is connected.
- 5 Check the key switch of machine is connected.

## Chapter 4: Program setting

1. Press power button to start the machine, the screen will show “Welcome”.

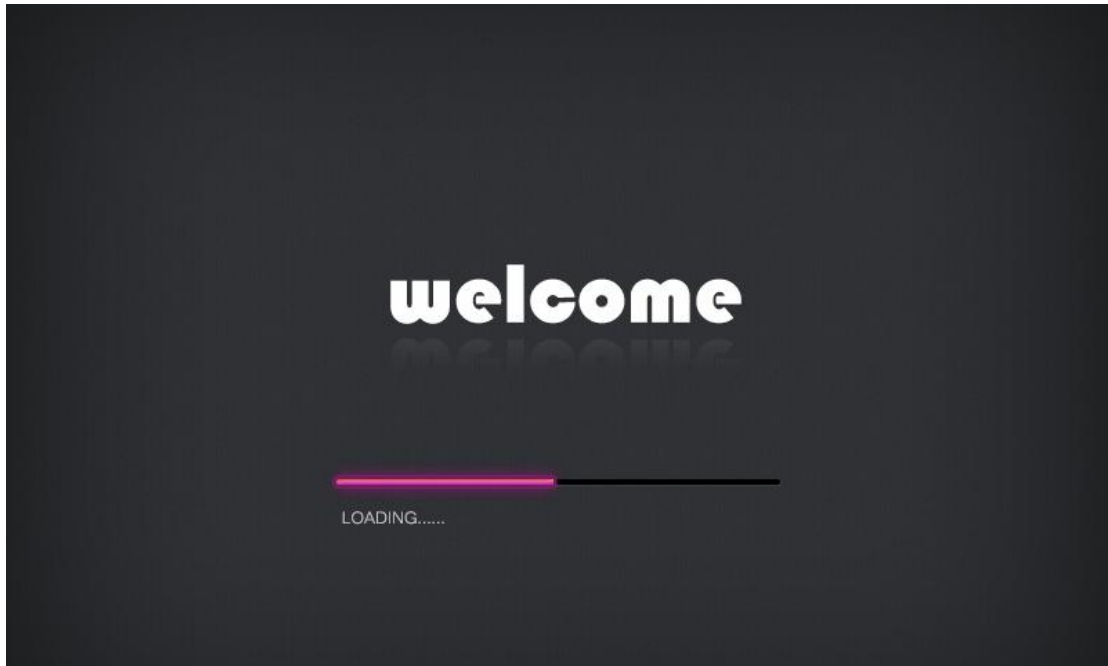


- 2 Press welcome to enter the main menu.

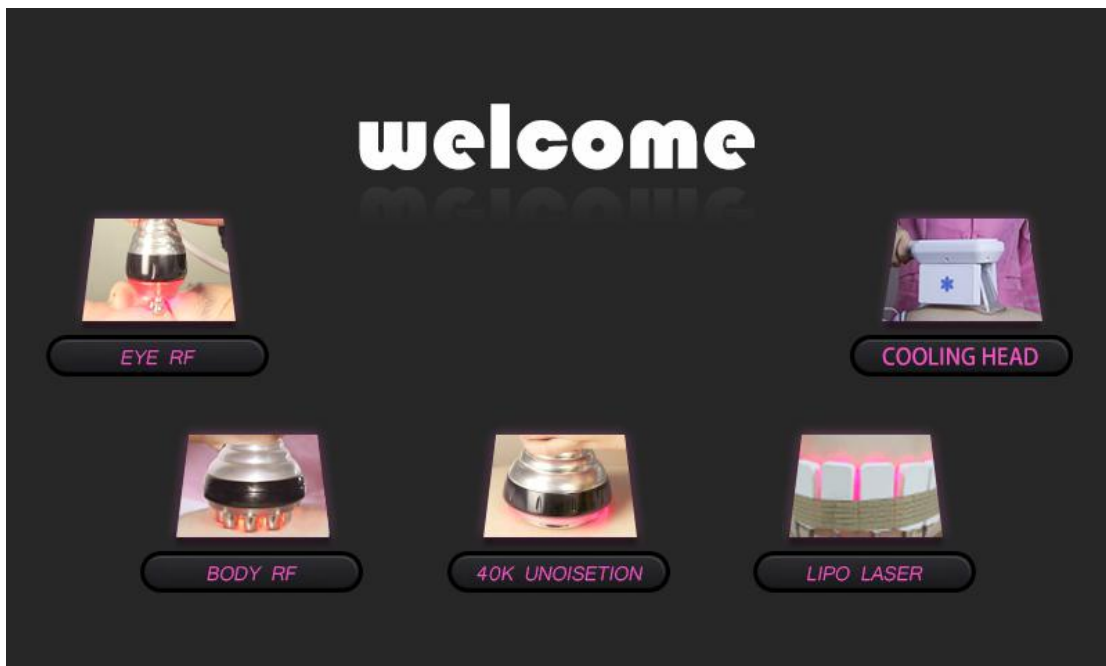

3. Choose the COOLING HEAD, and set parameters as below:



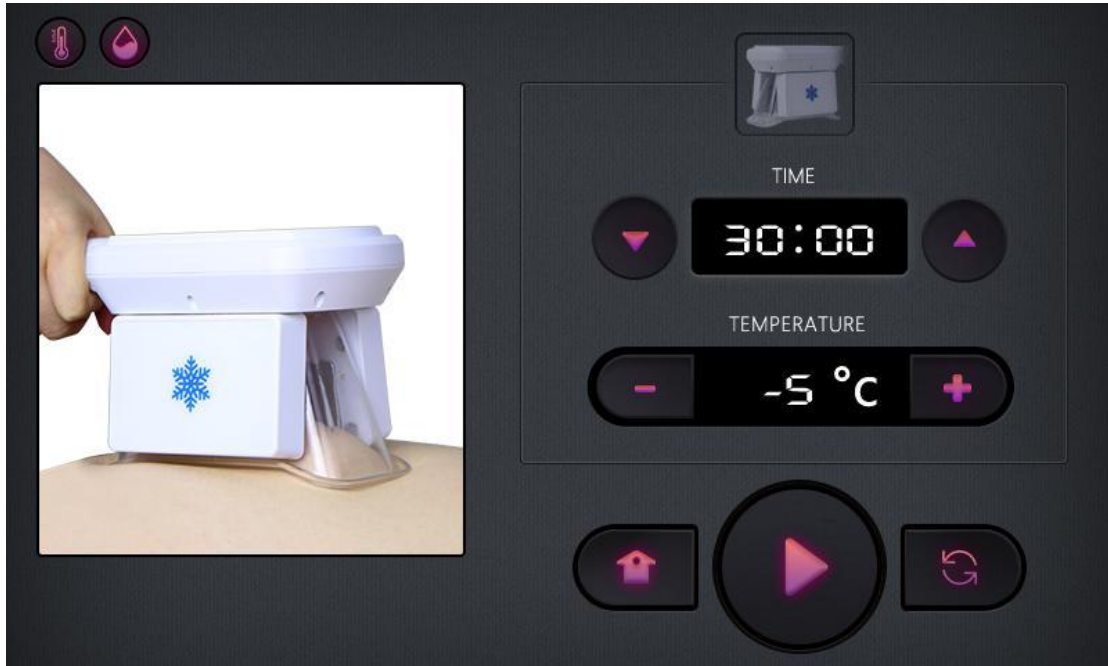

Notice:



it is the temperature controller. If the temperature is too high, the machine will stop itself and this icon will become gray



It is the water controller. If the water is not full, or there is something wrong with the water circulation, the machine will stop itself in 5 seconds to protect the machine, and this icon


will become gray 


A. Please click  to set working time


B. Please click  to set temperature. from high to low ,

gradually find the temperature suitable for the skin


C. Please click  to start or stop working

D. Please click  to return to home page


E. Please click  to clear the setting parameter

F. "Pressure" button  refers to suction intensity, which can be set as needed.

4 Choose the EYE RF, and set parameters as below:

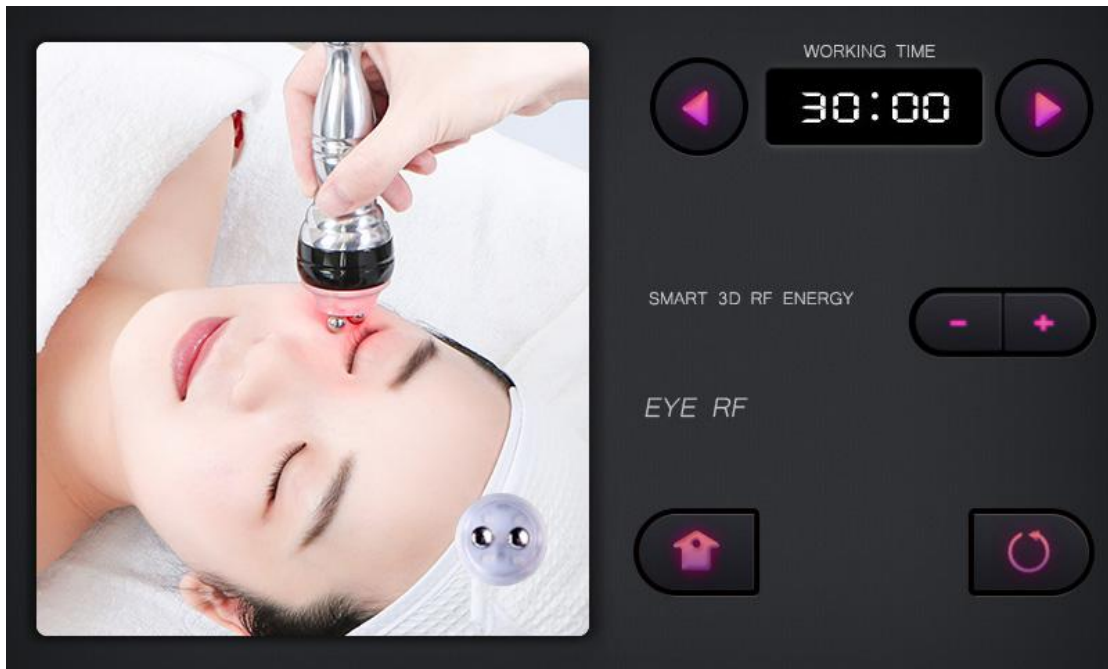

A. Please click  to set working time

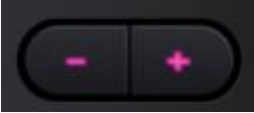
B. Please click  to set the energy. from low to high, and gradually find the temperature suitable for the skin




C. Please click  to choose mode.

NOR mode: the temperature directly reaches the set energy . lights up even if there is no contact with skin

PRO mode: the temperature gradually reaches the set energy . lights up when the head full contact with the skin

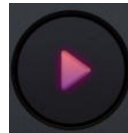


D. Please click  to start or stop working

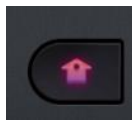


E. Please click  to return to home page




F. Please click  to clear the setting parameter




G. Please click  to switch the RF head of the red light


5.Choose the BODY RF, and set parameters as below:

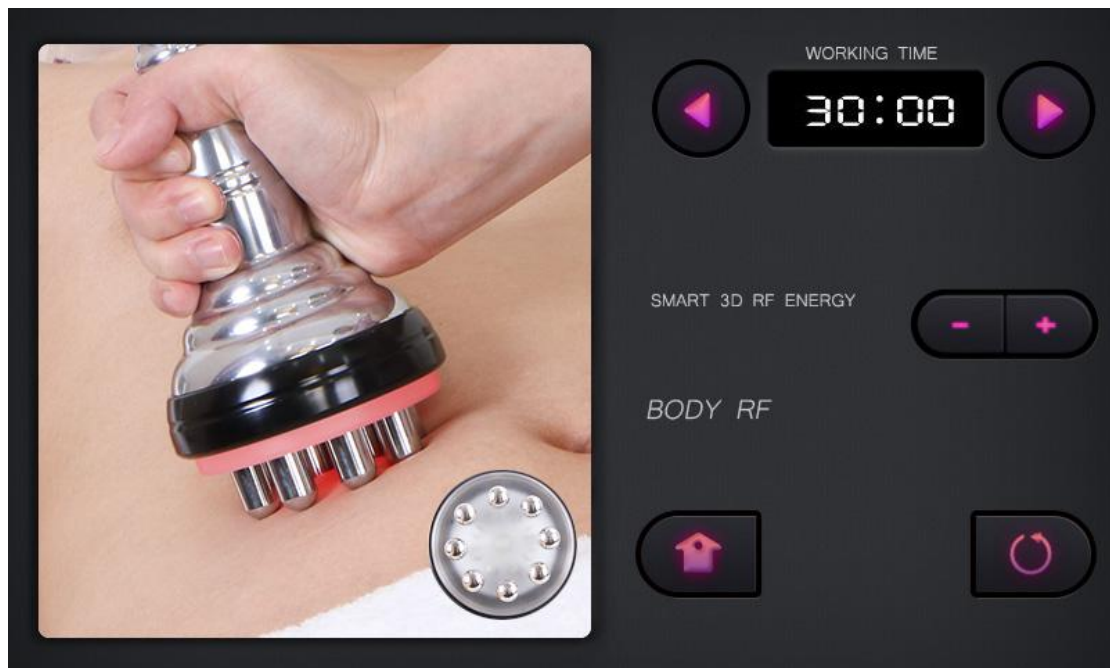




A. Please click  to set working time

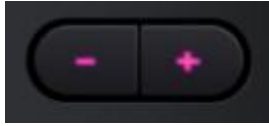

B. Please click  to set the energy, from low to high, and gradually find the temperature suitable for the skin



C. Please click  to choose mode.

NOR mode: the temperature directly reaches the set energy . lights up even if there is no contact with skin

PRO mode: the temperature gradually reaches the set energy . lights up when the head full contact with the skin




D. Please click  to start or stop working




E. Please click  to return to home page





F. Please click  to clear the setting parameter


6.Choose the 40K UNOISETION, set parameters as below:

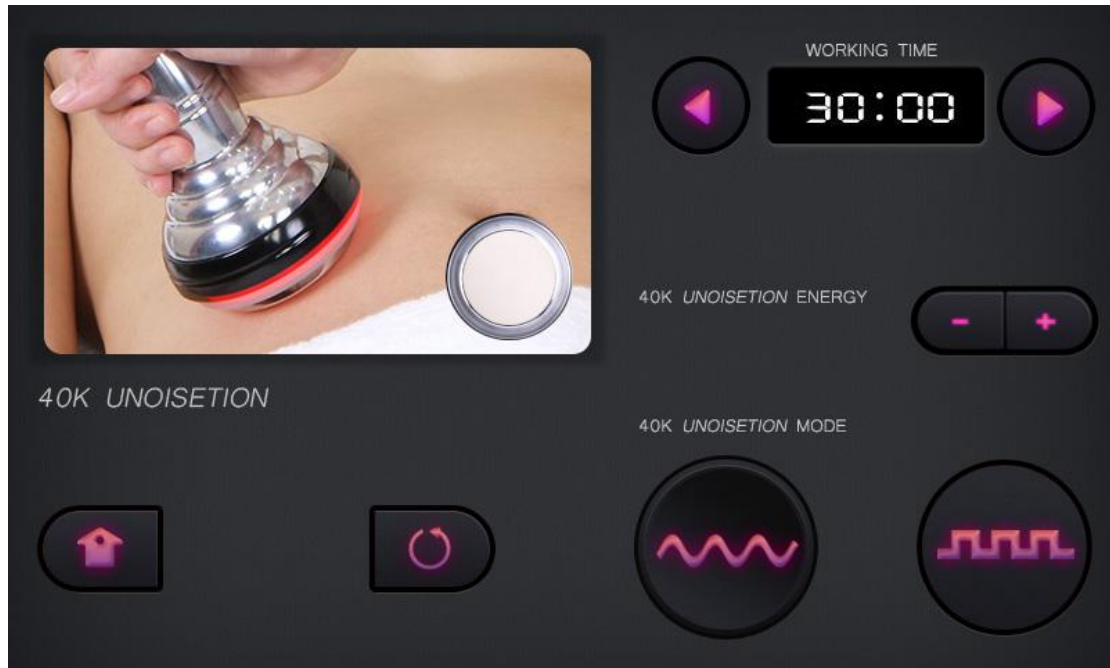


A. Please click  to set working time

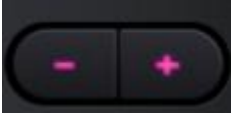B. Please click  to set the energy

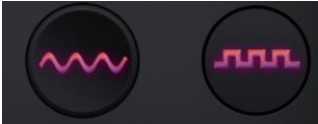C. Please click  to set mode

D. Please click  to start or stop working


E. Please click  to return to home page


F. Please click  to clear the setting parameter

7. Choose the LIPO LASER, set parameters as below:





A. Please click  to set working time

B. Please click  to set mode

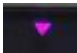
C. Move the icon  to the laser pad which you need to adjust the energy.

Please Click  to set the energy.

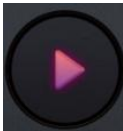D. Please click  to start or stop working

E. Please click  to return to home page



F. Please click  to clear the setting parameter

## Chapter 5: Treatment order

### I Fat freeze body slimming(abdomen treatment for sample)

- A. Disinfection: please use the cotton with medical alcohol to disinfect the treatment area
- B. Apply 2 pieces of antifreeze membranes on the treatment area
- C. Set the working time and temperature of cooling handle, the temperature range is 15~ -10°C, click to start
- D. Adjust the vacuum pressure.Put the handle on the treatment area and click start button, handle is working.
- E. Turn off vacuum, put the handle on the holder
- F. Take away the antifreeze members
- G. Recommendation

Make sure set the temperature from high to low upon the feeling of customer , -2 to 5 °C recommended

Treatment recommendation: 20-30minutes/session, 1 session/month, 3 session/course of treatment

### II Unoisetion2.0 cavitation+RF+Lipo laser body slimming

(abdomen treatment for sample)

- A. Disinfection: please use the cotton with medical alcohol to disinfect the treatment area
- B. Apply slimming cream on the treatment area
- C. Set the working energy and time of body RF probe, use the probe to massage the specific part of skin circularly after applying the slimming cream to promote activity of fat cells

Treatment time: 5 minutes recommended

D.Set the working energy, time of cavitation probe and mode, two modes can be chosen. Then use the probe to massage the specific part of skin circularly to cause effective consumption of calories and moisture in fat cells and reduce the size of fat cells. What's more, make fat cells be exploded instantaneously and disintegrate triglyceride into glycerol and free fatty acids, reduce the amount of fat cells and thereby achieve the effects of removing fat.

Treatment time: 5-10minutes recommended

D. Set the working energy and time of body RF probe, use the probe to massage the specific part of skin circularly to enhance blood circulation, accelerate the decomposition of fat, take away wastes such as glycerol and free fatty acids, enhance the effect of removing cellulite

Treatment time: 5 minutes recommended

E. Make the lipo laser pads on the specific part of skin by pairs of band, and set the time,output energy of lipo laser, to help body take away wastes more fastly, at the same time,it can help to



make fat cells reduce too

### **Skin lifting**

#### **EYE RF**

- A. Massage the head to relax body and face
- B. Please use the cotton with medical alcohol to disinfect hands
- C. Skin cleaning (less than 3minutes)
- D. Use the probe to massage the specific part of skin circularly after applying the anti-aging essence
- E. Set the working energy and time of face RF probe, use the probe to massage the specific part of skin according to the facial muscle distribution by beautician's skin

Treatment time: 5-10 minutes recommended

- E.Do moisturizing mask on face with LED photon therapy

Treatment time: 5-10 minutes

- F. Do moisture plaster mask

Treatment time:5-10 minutes


- G.Finish cleaning

- H.Do daily skin care

Treatment recommendation: 30-40minutes/session, 3sessions/week, 10 sessions/course of treatment, 2-3 courses of treatment will be more better

## **Chapter 6: Treatment skills**

### **Cryolipolysis**

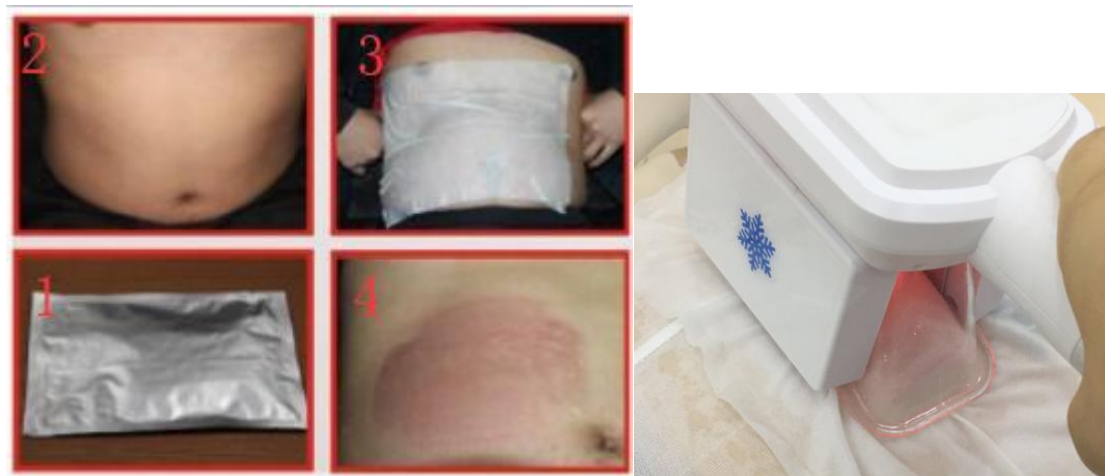

### **40K UNOISETION, BODY RF, EYE RF**

#### **THIGH**

Treated time: 30 minutes





1. From down to upper, pushing to the groin to dredge the lymph.

2. From down to upper, circling push by anticlockwise can help decomposing



fatness.



3. Also can push by come-and-go to decompose fatness.

4. Pull from the knee and muscle texture to upper, can improve the curve.

Treated time: 20 to 30 minutes

1. Along the lymph direction, pulling to the waist.



2. From upper to down, by anticlockwise gesture, pull come-and-go.



3. Along the muscle of arm, pull up to lift and tight the muscle.

## BACK

Treated time: 20 to 30 minutes



1. First by come-and-go, pushing the bladder nerve 2 to 3 times.




2. Circling by anticlockwise to stimulate the underarm lymph node.



3. On the back, by anticlockwise circling, can help to decompose the fatness.

4. Using lymph drainage gesture, pull the toxin to the lymph node, can help improving the back

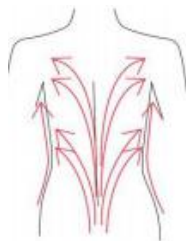

curve.

## ARM

Treated time: 30 minutes

### BACK SIDE

1. From the inner elbow to the armpit, doing lymph drainage.
2. Circling to stimulate the lymph node.
3. By anticlockwise or come-and-go pushing, decomposing the fatness.
4. From the elbow to the oxtter, doing tightening gesture.





## Facade

1. From the elbow to the oxtter, doing tightening gesture.

## Abdomen

Treated time: 20 to 30 minutes

1. From small to large, by clockwise direction, circling around the navel, can help peristalsis of the large intestine.
2. From small to large, by anticlockwise direction, circling push down.
3. From the belly and muscle texture, pull to the groin.
4. Through lymph drainage gesture, taking the toxin to the place of groin.

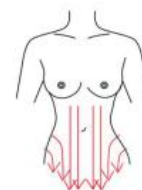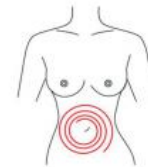

### Attention:

1. Pregnant women or women during in menses.

## Chapter 7: Attention for cryolipolysis

1. The treatment time at each body site is 20-30 minutes (depending on the thickness of the fat and health status).
2. Eliminate body fat part by part, i.e. it is required to complete a few courses to eliminate the fat of the entire waist.
3. After a treatment about 30 days, take another cryolipolysis treatment at the same location, which can reduce the excess fat thickness. Such a secondary treatment can consolidate the fat dissolving effect, but multiple treatments are not recommended.
4. Cryolipolysis cool-sculpting treatment is designed for those physically healthy people who have local fat deposits, but not suitable for the one with overall obesity. It can target at the fat at a specific part, such as waist, abdomen, butt, etc. The treatment is not a panacea, it can not treat obesity. It is not a substitute for liposuction or other invasive treatments. Substantial weight loss should be achieved through a balanced diet and regular exercise, cryolipolysis cool-sculpting treatment can only eliminate the extra fat.
5. Qualified antifreeze fluid or Antifreeze membranes are necessary to avoid frostbite during the treatment process when the temperature is reduced to 0 to -40 degrees. [Consumables, purchased by the user], Please apply at least 2pcs membranes on the treatment area
- 6 Use appropriate and constant suction intensity, the suction intensity is adjustable according to the patient's individual vulnerability. At the same time, the mode of suction and massage can be used.
- 7 After 25-30 minutes of treatment, press the pause button to suspend the suction, check and add antifreeze fluid.
- 8 Considerations after treatment

After the treatment, skin injury, minor frostbite, or pigment may occur, pay attention to take good care of it.

## Chapter 8: Contraindications

- 1 The patient cannot receive the treatment if he/she has following medical electronic equipment:
  - Implanted cardiac pacemakers or other implanted electronic devices;
  - Life-support artificial heart-lung machine or other equipment;
  - Portable ECG measurement device;
- 2 The following patients shall have the professional physician's consent before using the machine

- Acute disease patient; - Cancer patient; - Infectious disease patient; - Pregnant women;

- Heart disease patient; - Fever patients;

- Patients receiving another treatment or those with physical anomalies;

3 Careful attention should be paid to following patients when using this machine:

- Infants and young children; - People with sensory or physical impairment;

#### 4. Contraindications

Liver problems; Raynaud's disease; severe cold urticaria; cryoglobulinemia; cold fiber hypoproteinemia; severe frostbite; severe diabetes; as well as the elderly, children, the frail, and the one extremely vulnerable to cold

Product name :		Purchase date:		deadline:	
Seller:		Add.:		Tel.:	+86 18948755929
Customer:		Add.:		Tel.:	
Repair date		Reasonof the fault		Repair content	
Maintain list					