

WL-UU551

User manual



Catalog

Warning.....	2
Parts instruction	3
Accessories	4

Basal operation	5
40kcavitation Explode procedure.....	5
Caution for using 40k.....	6
Calibration for 40k.....	6
Sixpolar RF for body	6
Tripolar RF for face.....	7
Bipolar RF for eye.....	8
Caution for using RF	8
Forbidden group.....	8
Technical parameters.....	9
Packing listing	9
FAQ.....	10

Warning

Dear user ,

Thanks for using our company's products,in order to fully make use of the products , we sincerely suggestion follows:

- 1.Please read the instruction in detail and keep it well for further reading and consulting.
- 2.Please follow the instruction to fit and operate the machine correctly.

3. Please do not remove and change any accessories of the machine.
4. It forbids others to open and remove the mainframe except the authorized people in the company .
5. It forbids others to change and maintain the machine except the professional people.
6. Please contact us if there is any trouble during operation.

Parts Introduction





connecting sockets

MYCHWAY® 专利产品
仿冒必究
PATENT REGISTRATION PIRATED RESERVED

Accessories



40k cavitation head




Sixpolar RF + Laser



Tripolar RF Head



Bipolar RF

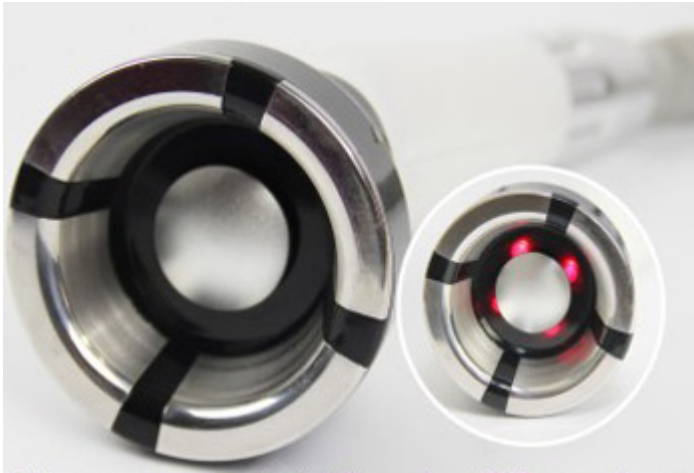

Vacuum + Multipolar RF

Basal Operation


1. Put all the accessories well and connect them to the host .

- ① Insert the Sixpolar laser RF work head to (A)
- ② Insert 40k cavitation work head to (B)
- ③ Insert Tripolar RF work head to (C)
- ④ Insert Bipolar RF work head to (D)

2. Insert power supply to (F), Turn on general power switch (E)



Interface LCD screen

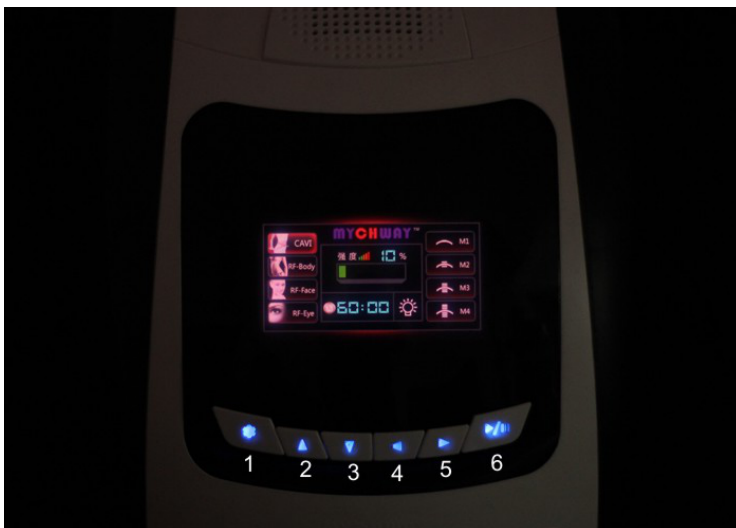

40k cavitation Explode Procedure

Apply the cavitation theory ,use the 40Khz Ultrasound to make intense vibration , which causes a inner explosion in the fat issues , thus get rid of lipocyte, it can effectively break down stubborn and deep fat, decrease lipocyte quantity and its size.

The diathermal of the low frequency ultrasound can promote lipocyte and tissues liquid flow, which help drainage of fatty acid, waste and toxin through lymph system,

1. Clean the parts and applied gel on fleshly parts (e.g abdomen,hip, thigh), the amount depend on moving flexibility of the cavitation explode,
 - (1),Click the <1>, to choose the mode 40k cavitation ., time,mode and energy will be displayed
 - (2),Energy adjust .Click <4>,<5> for adjusting energy, <4>for decreased, <5> for increased.
 - (3), Time adjust .Click <2>,<3> for adjusting energy, <3>for decreased, <2> for increased.
 - (4), Mode Transfer,Click <6> first , make it to the interface for mode transfer , then click <1> to choose the mode which you like. When we do every mode, we can click <6> to let the laser light turn on and off.
- 2,Beautician hold the cavitation explode work head to move slowly on skin, in circle or beeline repeatedly , and the other hand push the fat towards the work head.
- 3.Do not use on back, avoid using on bones.
- 4.Do not use on uterus part when a woman do abdomen treatment during menstruate should be avoid .
- 5.The machine will pause automatically if the setting time is over , Please click any place the screen according to the hint.
- 6.Clean the leftovers with the hot towel ,use warm towel to clean the cavitation explode work head.

Caution for using 40k

- 1.Applied enough gel on skin need treated
- 2.Do not stay unmoved on any part of your body , do not stay on bones.
- 3.The total treatment time for daily should be within 30minutes.
- 4.There will be ultrasonic only when you start the working mode of 40k.
- 5.Do not allowed to be placed under the state of no-load to avoid the damage caused by high heat.
- 6.40k work head can be used on fresh part of the body, Do not use for facial , If the skin is sensitive or becomes unhealthy

Calibration for 40k

1. We can click the <1>and <6> to calibrate the frequency of 40k cavitation head, so that get the best result and comfortable ways to do the treatment.

Sixpolar RF for body

- 1.Click <1>, Adjust the mode to Sixpolar RF,.
 - (2),Energy adjust .Click <4>,<5> for adjusting energy, <4>for decreased, <5> for increased.
 - (3), Time adjust .Click <2>,<3> for adjusting energy, <3>for decreased, <2> for increased.
 - (4), Mode Transfer,Click <6> first , make it to the interface for mode transfer , then click <1>

to choose the mode which you like. When we do every mode, we can click <6> to let the laser light turn on and off.

2. Set works time, 20 minutes for each parts

3.The system default laser is start state , click <6>, the laser function is closed.

Note:Do not let the laser light straight shine upon the eyes in case of eyes hurt

4.Applied gel on skin needed nursing , the electrode probe move slowly on skin, if there is too much fat , you can increase the output intensity.

5.The machine will pause automatically if the setting time is over .

6. Clean the leftover with hot towel,use warm towel to clean the work head.

Tripolar RF head for face

1.Click <1>, Adjust the mode to Tripolar RF,.

(2),Energy adjust .Click <4>,<5> for adjusting energy, <4>for decreased, <5> for increased.

(3), Time adjust .Click <2>,<3> for adjusting energy, <3>for decreased, <2> for increased.

(4), Mode Transfer,Click <6> first , make it to the interface for mode transfer , then click <1> to choose the mode which you like. When we do every mode, we can click <6> to let the laser light turn on and off.

4. Set works time, 20 minutes for each parts

5.The system default laser is start state , click <6>, the laser function is closed.

Note:Do not let the laser light straight shine upon the eyes in case of eyes hurt

6.Applied gel on skin needed nursing , the electrode probe move slowly on skin, if there is too much fat , you can increase the output intensity.

7.The machine will pause automatically if the setting time is over .

8. Clean the leftover with hot towel,use warm towel to clean the work head.

Bipoar RF head for eye

1.(1).Click <1>, Adjust the mode to Bipolar RF,.

(2),Energy adjust .Click <4>,<5> for adjusting energy, <4>for decreased, <5> for increased.

(3), Time adjust .Click <2>,<3> for adjusting energy, <3>for decreased, <2> for increased.

(4), Mode Transfer,Click <6> first , make it to the interface for mode transfer , then click <1> to choose the mode which you like. When we do every mode.

2.Set works time, 20 minutes for each parts

3.The system default laser is start state , click <6>, the laser function is closed.

Note:Do not let the laser light straight shine upon the eyes in case of eyes hurt

4.Applied gel on skin needed nursing , the electrode probe move slowly on skin, if there is too much fat , you can increase the output intensity.

5.The machine will pause automatically if the setting time is over .

6. Clean the leftover with hot towel,use warm towel to clean the work head.

Caution for RF

1. Beautician should wear insulated gloves during operation.
2. Do not wear any aglet, both for the beautician and the users
3. The electrode probe must be fully in contact with skin.
4. Do not make the RF electrode stay on any part without moving, Keep it moving on skin slowly.
5. Skin around eyes and forehead is very thin, the RF output intensity should be relatively weakened.
6. Do not contact skin with unipolar when using of bipolar, tripolar, hexapolar work head otherwise it will produce pulse discharge, gel paste. Ect. Products like gel maybe evaporated into a mist.
7. Do not let the laser light straightly shine upon the eyes in case of eyes hurt.

Vacuum Multipolar RF

1. (1). Click <1>, Adjust the mode to Vacuum.
- (2). Energy adjust. Click <4>, <5> for adjusting energy, <4> for decreased, <5> for increased.
- (3). Time adjust. Click <2>, <3> for adjusting energy, <3> for decreased, <2> for increased.
- (4). Mode Transfer, Click <6> first, make it to the interface for mode transfer, then click <1> to choose the mode which you like. When we do every mode.
2. Set work time, 20 minutes for each parts
3. The system default laser is start state, click <6>, the laser function is closed.
Note: Do not let the laser light straight shine upon the eyes in case of eyes hurt
4. Applied gel on skin needed nursing, the electrode probe move slowly on skin, if there is too much fat, you can increase the output intensity.
5. The machine will pause automatically if the setting time is over.
6. Clean the leftover with hot towel, use warm towel to clean the work head.

Forbidden Group

1. People with heart disease or high blood pressure, or configured cardiac peacemaker.
2. Patients with acute inflammation, asthma, deep vein thrombosis, thyroncus, cancer.
3. People with hemorrhagic disease, trauma or who is bleeding.
4. Pregnant women
5. Medical Plastic parts, or parts with artificial in filling inside.
6. People with metal inside the body (include metallic tooth)
7. Patient with skin inflammation or with edema.
8. People with an abnormal immune system.
9. Numb or insensitive to heat.

Technical Specification

Screen:5.4" Color TFT Display

Rated input voltage:AC220V-240V 50HZ /AC110V-120V 60Hz

Rated input power:350VA

Cavitation frequency : 40KHz

Cavitation Power:up to 60W

Cavitation tip:60mm diameter

RF for body

RF frequency:3MHz

RF power:up to 150W

RF type:six polar

RF tip:65mm diameter

RF for face

RF frequency:1MHz

RF power:up to 100W

RF type:tripolar

RF tip:25mm diameter

RF for eye

RF frequency:1MHz

RF power:up to 60W

RF type:bipolar

RF tip:25mm diameter

Lipolysis laser:635nm

Lipolysis power:100mw

Photon

Light:Red LED light

Wavelength:640nm+-5nm

Power:15W

Packing Include

Main machine x1 unit

40KHz cavitation head x 1unit

Sixpolar RF head for body x1unit

Tripolar RF head for face x 1unit

Bipolar RF head for eyes x 1unit

Power supply wire x 1unit

Holder x 1unit

Key switch x 1 unit

FAQ:

1.How quickly will I see the results of treatment?

Most clients experience few cm of reduction in circumference after a single session, with increasing results after each visit. The results may vary with different tissue structure, treatment area, age, metabolism, medications, and changes in hormones. Proper diet and increased physical activity will certainly improve and help to maintain the results.

Between 8 and 12 sessions, you will be able to appreciate the results from the third session on.

2.Is repeat treatment required?

Generally, one treatment per a single area is usually adequate to achieve the required results. And we suggest that you had better have 8 and 12 sessions.

Normally it is a Long lasting treatment.

However, maintenance of this fat loss is dependent on a healthy diet and exercise.

3.How does it feel?

Most clients consider the procedure painless and comfortable. There may be, however, a slight discomfort due to the specific noise experienced, but it poses no harm and disappears as soon as you are not in contact with the applicator. You may also experience a little warmth during the treatment.

It is very normal if you experience with this above problems in which is doing the treatment of Real 40Khz cavitation.