

# User Manual

SR-002



# Preface

Dear Users,

We're pleased to present to you our latest beauty machine SR-002, which is a skin-cleansing and acne removal product based on high-frequency electromagnetic wave technology. It is a professional beauty machine, which needs to be operated by professionals with proper training. Any improper use will result in adverse outcomes. Therefore, we advise people to read this manual thoroughly and follow the instructions strictly before operation.

We believe our quality product will bring in good returns on your end. Also, you can count on our worry-free after-sales service.

Thank you!

# Table of Contents

Preface	
Part I	
Introduction	
Advantages	
Working Principle	
Skin Cleansing&Acne Removal	
Head Treatment	
Part II	
1. Detailed Operation	
2. Technical Specifications	
3. Dos and Don'ts	
4. Troubleshooting & Solutions	
5. FAQs	
6. Packing List	
7. Simplified Treatment Steps	
8. Operational Diagrams	

# Part I

## Introduction


The beauty machines based on the high-frequency electromagnetic wave technology are widely used and very popular among people, which can effectively replace some traditional acne removal treatments in beauty salons. Moreover, it is convenient and easy to operate. It solves acne-related problems for beauty-seekers. This equipment requires no injection, no medication, and no surgeries, for which it has no side effects. It is used externally during the whole process and has an immediate effect.

## Advantages

1. Working tips with multiple shapes can be chosen for use according to different needs
2. Ozone can be produced under the influence of high-frequency electromagnetic wave, which has the ability of relieving inflammation and killing bacteria.
3. It is suitable for acne-prone skin, which can relieve acne and repair acne scarring at the same time.
4. Comb-shaped working tip can help reduce dandruff and hair texture.
5. It is painless and non-invasive treatment process.
6. Once the treatment is over, customer can go back home and there is no recovery period, which brings no inconvenience to customer's life and work.
7. No consumption, low cost and quick returns.
8. Unevenness, bleeding, and swelling will not appear after treatment.

## 1. Working Principles

High-frequency electromagnetic wave can promote the vitality of basal layer of skin, improve skin texture, increase skin's resistance to bacteria. Thanks to the production of ozone, it can help kill bacteria. It can limit sebum secretion, help with the healing process of infected skin. It can also help open up pores to allow impurities to be removed easily. It can promote lymphatic and blood circulation, promote skin's texture and regulate the balance between moisture and grease. The Comb-shaped working tip can promote qi and circulation on head, and prevent hair loss and gray hair.

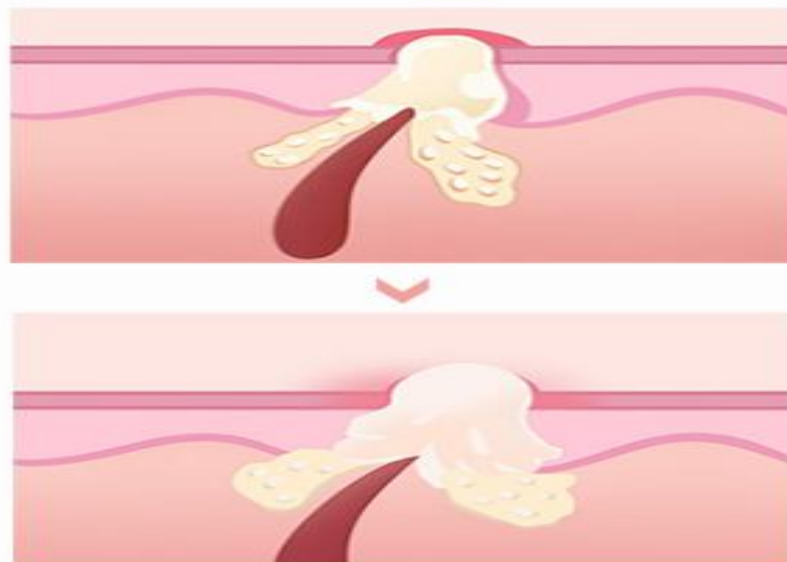

## 2. Skin Cleansing & Acne Removal

### Effects

1. Relieve inflammation and kill bacteria for deep skin layer.
2. Relieve acne, fade spots, whiten skin and rejuvenate skin.
3. Relieve skin's roughness and large pores.
4. Repair acne-damaged skin.
5. Help with removal of mites.
6. Promote circulation and increase skin's metabolism.
7. Brighten skin, minimize pore size.

### Applicable Range

1. Those who have active grease secretion and large pores.
2. Those who have acne-prone skin.
3. Those who have large-area acne.
4. Those who suffer from acne-related facial problems.
5. Those who suffer from blocked pores and mite-infested skin.

### Inapplicable Range

1. Those who have severely sensitive skin or are having allergic reactions.
2. Those who have prostheses or have large-size wounds on face.
3. Pregnant women or patients with severe heart diseases, severe hypertension, hyperglycemia, hyperlipemia, major diseases and so on.
4. Those who have skin diseases or contagious diseases.
5. Those who have just taken skin booster, and micro-needle or who have wounds on skin.
6. Those who have unrealistic expectations about the effects.

### Do's and Don'ts after Treatment

1. Drink more water and use sunscreen.
2. Do not use irritating products such as AHA, scrubbing and exfoliating products within one week.
3. Do not stay up late, avoid drinking alcohol and smoking cigarettes and eat more fruits, vegetables and foods with less grease.
4. Apply face mask three days in a row.
5. Wash face with hot water within 3 days.
6. Binge eating is not allowed and do not eat spicy foods.

## 3. Head Treatment

### Effects

1. Promote qi and circulation on head.
2. Promote hair growth.
3. Help reduce dandruff, prevent hair loss and oily hair.
4. Relieve the condition of hair dryness and proneness to break.
5. Improve skin texture and prevent gray hair.
6. Relieve insomnia, too many dreams in sleep, dizziness and headache.

## Applicable Range

1. Those who often suffer from dizziness, and the condition of there being too many dreams in sleep.
2. Those who often have itchiness on their scalps.
3. Those who suffer from hair loss and whose hair is prone to break.
4. Those who have poor hair texture and split ends.
5. Those who have dry and oily hair.
6. Those who have very thin hair.

## Inapplicable Range

1. Those who have super sensitive skin.
2. Those who suffer from ulceration and infection on scalp and who suffer from contagious skin diseases.
3. Pregnant women or patients with severe heart diseases, severe hypertension, hyperglycemia, hyperlipemia, major diseases and so on.
4. Those who are of advanced years.
5. Those who have just taken operations on scalp and are on the road to recovery.
6. Those who have unrealistic expectations about the effects.

## Do's and Don'ts after Treatment

1. Do not wash hair immediately after head treatment.
2. Do not wash face or hair within 6 hours after treatment.
3. Do not expose yourself directly to the wind from electric fan or air-conditioner.

## Part II

### 1. Detailed Operation

Choose proper working tip  and install it on the handpiece . As the red arrow shows , you can see the power switch and button for adjusting energy level there. Adjust energy level based on individual endurance.

### The functions of different tubes

**Comb-shaped tube:** used on head. Suitable for those who have poor hair texture and suffer from insomnia and the condition of there being too much dreams in

sleep.

**Crooked-shape tube:** suitable for nose treatment to relieve acne and kill bacteria

**Spoon-shaped tube:** suitable for face treatment to relieve inflammation and acne on small area

**Mushroom-shaped tube:** suitable for face treatment to relieve inflammation and acne on large area

## 2. Technical Specifications

110v~240v

## 3. Dos and Don'ts

1. The plug with grounding pin must be used before the instrument is used and the power socket of the instrument is properly grounded.
2. Ensure that the voltage of the instrument is suitable. If the voltage of the local power supply is not stable, we suggest that the user add a voltage regulator with matching power between the mains and the instrument.
3. In order to ensure the treatment effect and normal service life of the instrument, please use the specified parts provided or suggested by the original manufacturer.
4. The instrument should not be placed in damp places or near water sources, and should not be directly exposed to sunlight.
5. Do not place the instrument near a strong heat source, as this may affect the life and normal use of the instrument.
6. Before treatment, please remove all metal objects from the body first, so as to avoid unexpected conditions and influence the curative effect.
7. Please do not use the instrument on eyes, thyroid, parathyroid, testicles, abdomen of pregnant women, and pacemakers.
8. Patients who are suffering from illness should consult a doctor and get permission from the doctor before using.
9. Please turn off the power switch of the device, and ensure that the total power supply of the device is turned off before the staff leaves after daily use, so as to ensure the safety of the electrical products.
10. Use device or train device operators in strict accordance with instructions in the user manual.
11. Set the energy at the lowest level at the beginning of operation and increase the level slowly after adapting to it.
12. Do not use the instrument in eyes, ears and the areas where there is skin diseases or plastic surgery is taken.
13. Clean the instrument with normal saline after operation to ensure its cleanliness and hygiene and prolong its service life.

## 4. Troubleshooting & Solutions

### 1. The device cannot be started?

- A. Make sure the power cord is connected to the socket with normal power supply.
- B. The power cord is broken.

### 2. The handpiece has no energy output?

A: Please check whether the power switch is pressed.

B: The handpiece is broken, and you need to replace it with another one.

C: Please check whether the working tip is properly connected to the handpiece.

## 5. FAQs

### 1. Is there any dependence on the treatment?

A: No. This treatment only helps skin absorb nutrients and solves skin problems such as acne, pimples and roughness. If you don't take the treatment in later stage, skin will just return to its normal secretion and ageing process. There will be no dependence.

### 2. Can I wear makeup after the treatment?

A: it is recommended that you should not wear makeup immediately. This device can help the healing process of acne-damaged skin. The post-treatment makeup will be a hindrance to this process and even make skin get reinfected.

### 3. How long can I see results?

A: It depends on individual situation. The reasons why a person suffers from acne are different from each other's. Acne can occur due to both endogenous and exogenous factors. The length of treating acne is determined by the severity of acne, the self-restoration, regular management, and daily skincare habits. It is a gradual process with the treatment results having an accumulative effect.

### 4. Whether I will have acne again after the treatment?

A: there are many factors that contribute to your suffering from acne including the aspects of lifestyle, diet and skincare. In order to prevent acne, you can do something. Clear the excess grease out of follicles to keep them unblocked. Regulate the balance between moisture and grease.

## 6. Packing List

Purple Color High Frequency Hand Piece X 1

Mushroom-shaped Tube for large facial areas X 1

Spoon-shaped Tube for small facial areas X 1

Crooked-shape Tube for spot areas X 1

Comb-shaped Tube for Hair X 1

Instructions Sheet

## 7. Simplified Treatment Steps

### Skin-cleansing and acne removal




Remove makeup--clean face--apply blackhead-removing liquid and hot steam--pop pimples with a sanitized needle--machine treatment--apply face mask--apply water, essence, lotion/face cream, and sunscreen.

### Head treatment

Massage techniques---machine treatment



## 8. Operational Diagrams

Parameter Adjustment	Product	Techniques	Diagrams
<b>Skin Cleansing and Acne Removal: 45 min once a week</b>			
Energy level: rotate the lower part of the handpiece to adjust energy level.	Makeup-removal+Face cleanser+toner+face mask+SR-002	<ol style="list-style-type: none"> <li>1. Remove makeup, 2 minutes.</li> <li>2. Clean face with face cleanser, 3 minutes.</li> <li>3. Apply toner to the whole face, 1 minute.</li> <li>4. Apply blackhead removing liquid and hot steam( to the area where there are large pores or severe blockage. Toner can be applied to cheeks. If you have sensitive skin, apply toner and cold stem on the whole face instead of black head-removing liquid and hot steam)</li> <li>5. Use sanitized needle to pop pimples.</li> <li>6. <b>High-frequency tube</b> operation: choose the spoon-shaped tube, rotate the power switch to begin the treatment.</li> <li>7. Move the tube upwards in circular motion from chin in a line-by-line way, 3 times.</li> <li>8. Move the tube in circular motion from forehead to hairline, 3 times. (Repeat the operation on the area between the eyebrows 1-2 times. )</li> <li>9. Move the tube upwards from chin in a line-by-line way, 3 times.</li> <li>10. Move the tube from</li> </ol>	<p>Technique9</p>  <hr/> <p>Technique7</p>  <hr/> <p>Technique8</p> 

		<p>forehead to hairline, 3 times. (Repeat the operation on the area between the eyebrows 1-2 times. )</p> <p>11. Repeat the techniques on the other side.</p> <p>12. Choose crooked-shape tube, and then begin the treatment.</p> <p>13. Use the tube in a stamping manner on nasal fossa</p> <p>14. Choose the mushroom-shaped tube, use it a stamping manner on the whole face.</p> <p>15. Apply toner, 1 minutes.</p> <p>16. Apply face mask, 15 minutes.</p> <p>17. Remove mask and clean face, 2 minutes.</p> <p>18. Apply toner, essence, face cream and suncream, 3-5 minutes.</p> <p>19. The end.</p>	
<p><b>Treatment Suggestion</b></p> <p>One course of treatment consists of ten treatment sessions. After one treatment, pores will be visibly clean and acne condition will be relieved to some degree. After one course of treatment, skin tone will be gradually brightened up and pore size will be smaller with acne relived and acne scarring faded. Three courses of treatment, coupled with controlled diet, healthy lifestyle and good skincare habits, will prevent acne from recurring and fade acne scarring to the invisible level.</p>			
<p><b>Head treatment: 25 minutes twice a week</b></p>			
<p>Energy level: rotate the lower part of the handpiece to adjust energy level.</p>	<p>SR-002</p>	<p>1. Use thumbs to press the acupuncture points along Du Meridian alternately from hairline to the top of head, 2-3 times.</p> <p>2. Let head lie on its right side, use thumbs to press acupuncture points of the left side of head along Triple Energizer Meridian, 2-3 times.</p> <p>3. Do the same on the other side.</p>	

		<p>4. <b>High-frequency tube</b> operation: choose the comb-shaped tube, comb hair on the left side and then the right side, 3-5 times.</p> <p>5. The end.</p>	
--	--	--	--

**Treatment Suggestion**

One course of treatment consists of ten treatment sessions. After one treatment, dizziness, headache and migraine will be relieved and meanwhile circulation on head will be promoted. After one course of treatment, the results are more obvious. The possibility of the occurrence of dandruff, hair loss and oily hair will be drastically reduced. Three courses of treatment will help improve skin texture, and prevent gray hair and hair from breaking.