

# User Manul

The right to interpret this instructions belongs to the company



Product overview.....	01
Product Introduction .....	01
Product illustration.....	03
System Interface Diagram .....	04
Handle (A)Operation step .....	05
Handle (B) Operation step.....	07
V10 handle operation step.....	09
Operating methods .....	10
Use taboos.....	12
Attention.....	13
Troubleshooting .....	13
Instrument maintenance .....	14
Packing list .....	14

## Product overview

combined with HIFU/RF 40KHZ cavitation BIO microcurrent, vacuum, infrared multifunctional machine , it can get the effect of body slimming skin tightening . Please read this manual carefully before use and use the beauty equipment in strict accordance with the guidelines in this manual.

## Product Introduction

### ■ HIFU technology

Deep into the subcutaneous fat layer by turning the energy source into individual hot condensation points to destroy the unexpected fat cells, reduce unexpected fat cells, stimulate the contraction and thickening of collagen simultaneously, Stimulate Elastic fiber, unexpected fat cells decomposed into triglyceride, excreted through the body's metabolism to get the effect of body shaping .

### ■ RF Technology

With the technology of RF, skip the epidermis directly to the subcutaneous tissue, Resonating water molecules in the skin and subcutaneous tissue, accelerate heat generation by tissue friction , Elevated temperature promotes catabolism of subcutaneous fat, stimulating the proliferation and recombination of collagen and elastic fiber to achieve the effect of skin tightening , shaping and lifting.



## Product Introduction

### ■ V10 handle technologh

This handle is a combination of a varitey of multifunctional technologies. With the tech of 2\*40khz cavitaion, 2\*bio microcurrent, 2 RF , vacuum and infrared led . Firstly, using 40k cavitaion for fat explosion, Then use the high frequency rf act on water molecules in subcutaneous fat Precisely, Make bipolar water molecules produce rapid vibration, Heating the subcutaneous fat tissue by the heat generated by this high-speed rotary friction, Applying a special physical stimulus to the skin tissue to cause wrinkles on the compressed skin tissue, And subtle mechanical stimulation of deep skin tissue under the action of vacuum negative pressure, promote biochemical reactions in deep adipose tissue to accelerate fat consumption, rapid activation of lipase under specific magnetic and radio frequency to accelerate fat breakdown, decompose fat cells into free fatty acids and glycerol, to get the effect of decreasing volume of fat cells, improving cellulite tissue , body shapping.

## Product illustration











Power	160W	Voltage	220V
Operating Voltage	36V	Operation handle	3 pcs
Package size	52*46*28cm	Package weight	22kg



# System Interface Diagram

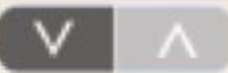

# Handle (A) Operation step

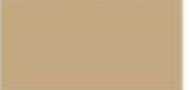





1. Clean the treatment area before using the machine. 
2. Connect the power cord and press the circular switch button at the top of the screen to turn on the machine. 
3. Click enter to enter the operation interface. 
4. Click handle (A) to enter the big handle operation interface,  Click handle (A) to enter the big handle operation interface,  for the first use the lowest energy 1 is suggested to use, then adjust the energy according to customer's tolerance, cause it is normal feeling hot or stinging when using the machine.
5. Click time to set operating time,  such as , the treatment time for waist is suggested about 30-45 minutes, it is depends on size of the operation area to set operating time.
6. After setting the operation time and operation energy , click the PAUSE,  There is a round button above the handle,  it is represent start or pause, When you press the button of the handle, you will hear the sound of "drop" then the machine starts work, when you need to pasue , press the button on the handle.



# System Interface Diagram

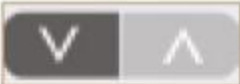


# Handle (B) Operation step


1. Click enter to enter the operation interface. 
2. Click handle (B) to enter the small handle operation interface,  The energy of HIFU is from 0-10, RF from 0-10, The up and down arrows is to increase or decrease energy,  for the first use the lowest energy 1 is suggested to use, then adjust the energy according to customer's tolerance, cause it is normal feeling hot when using the machine.
3. Click time to set operating time,  such as , the treatment time for waist is suggested about 30-45 minutes, it is depends on size of the treatment area to set operating time.
4. After setting the operation time and operation energy , click the  There is a round button above the handle,  it is represent start or pause, When you press the button of the handle, you will hear the sound of "drop" then the machine starts work, when you need to pause , press the button on the handle.



# System Interface Diagram

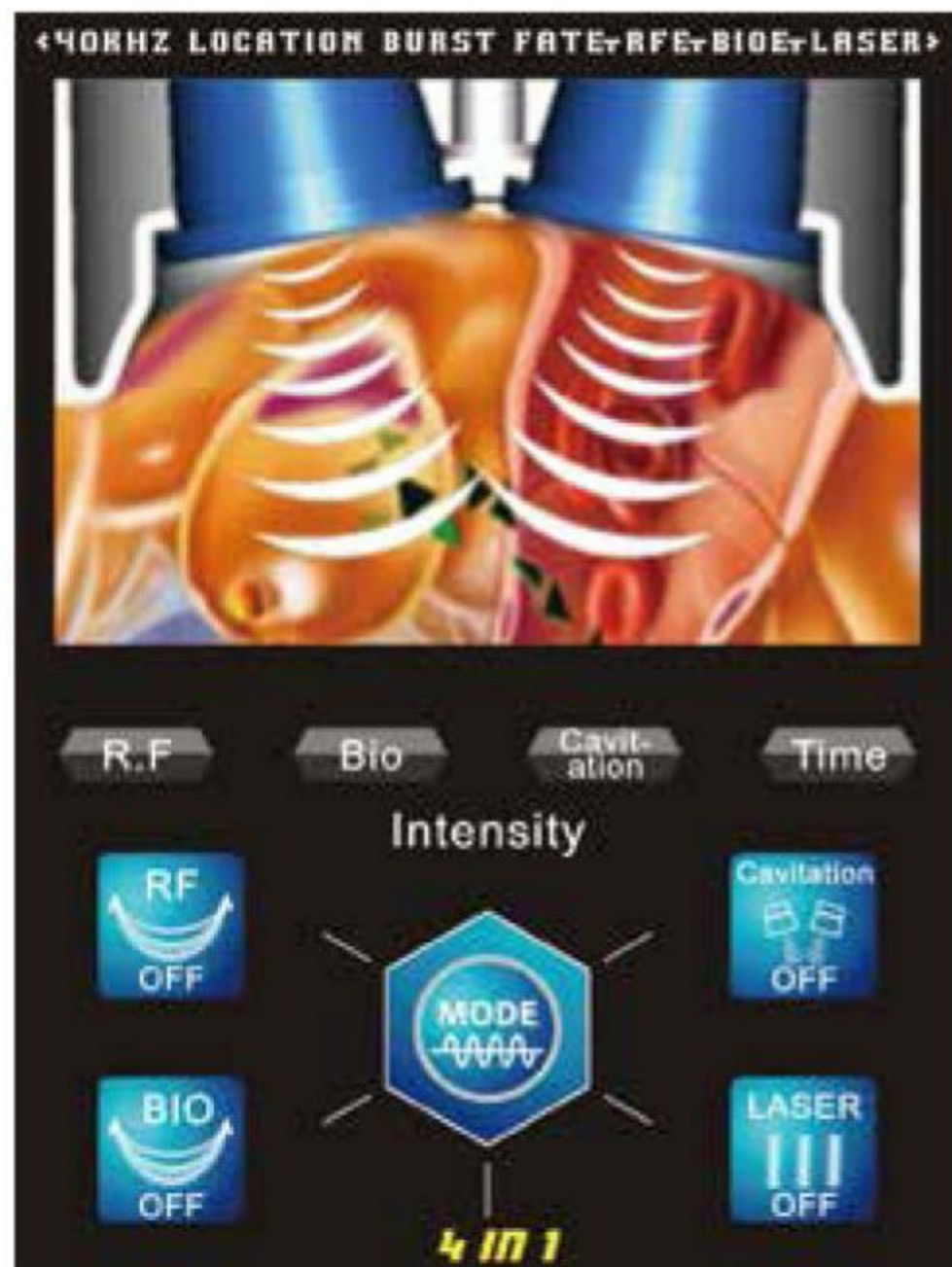

# V 10 handle operation step


1. Click "time" to set operation time,  it is suggested about 45-60 minutes, it is depends on size of the treatment area to set operating time.
2. Click energy setting, the energy of RF is from 1-10, bio is from 1-15, 40k cavitation is from 1-10, for the first use the lowest energy 1 is suggested to use, then adjust the energy according to customer's tolerance, It is normal feeling hot, Current sense or tinnitus when using the machine.
3. After setting the operation time and energy, Click the icon to turn red to represent output, and blue to stand by.
4. Clean the medium in the treatment area after the operation, then put on some moisturizing product.



# Operation method

1. Measure the size of the treatment area before using the instrument, Measure the marked area after the treatment.
2. Then apply the medium on the treatment area and the operating head.

## Waist operation

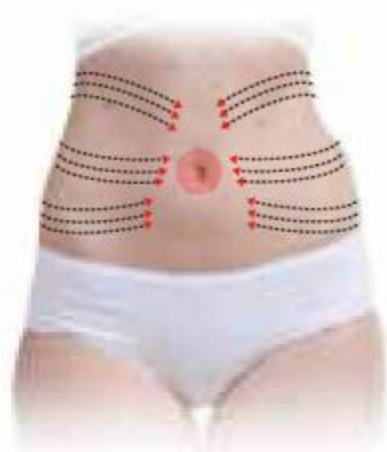

The operating head is first moved from the the upper waist to the navel around, then from the middle waist to the navel around, Lower waist moves to the navel, Finally moving from the navel to the groin , Do not move too large position during operation, it is recommended to move at a constant speed. Set the operating time according to the relaxation of individuals.

## Leg operation



The operation head moves from bottom to top in a circled motion, Move from the thigh armpit up to the thigh root , then from the top of the knee to the groin , Set the operating time according to the relaxation of individuals.

## Arm operation



The operation head is first pushed from bottom to top, lifting from the outside and moving to the vicinity of the armpit, Set the operating time according to the relaxation of individuals.

## Back operation



The operating head is operated from the underside of the cervical spine to the back waist, Going back and forth for the entire back.

## Hip operation



The operation head first operates from the bottom of the thigh and pushes from bottom to top;

## Backside waist



The operation head is placed on the upper hip bone, start operate from the position where the spine is about a fist.



## Use taboos

1. Women in menstrual period is not recommended to use.
2. Women in pregnancy or during the month after having a baby are not recommended to use.
3. Human who has metal, plastic, silicone such as implanted into body are not recommended to use.
4. Human who with hearing aid, Cardiac pulser such medical devices are not recommended to use.
5. Human who are hypertension, hyperlipidemia, hyperglycemia, or with metal problems are not recommended to use.
6. Human who are abnormal fat metabolism, including severe hyperlipidemia or severe fatty liver are not recommended to use.
7. Human with skin allergies or inflammation, skin infection are not recommended to use.
8. Human who use anticoagulant drug long-term continuously, such as: Non-steroidal anti-inflammatory drugs or contraceptive, long-term use of hormone drugs are not recommended to use.
9. Children under 12 years old or Human who has irrational expectation of the effect for the machine are not recommended to use.

## Attention

1. It is not recommended to take a shower with hot water immediately after 4 hours of operation;
2. After operation, it is recommended to maintain a replenishment of water up to 2000cc per day.
3. Keep at least 10-15 minutes of exercise per day, jog or do aerobic exercise is the best choice, keep a good mood.
4. Diet should be light, do not overeating, Keep nutrition matching.
5. Keep regular work and rest, do not over-drinking, eat too late at night, etc.

## Troubleshooting

Fault type	Troubleshooting method
The instrument did not respond after connecting the power	First check if the plug is powered, replace the plug or socket
Cannot enter the operation interface after turn on the machine	Shut down and restart, or loosen the screws on the screen

**Attention:** Instruments are precision electronic products, Using the instrument in a humid environment is forbidden: such as bathroom wet places; If it enter water accidentally, please stop using it and return back to the factory for repair!



# Instrument maintenance

1. The instrument must use a three-pronged plug with a grounding pin and ensure that the instrument's power outlet is well grounded;
2. Ensure that the power supply is stable and suitable. If the local power supply voltage is unstable, it is recommended that the user increase the power matching power supply.
3. Do not place the instrument in a humid place or near a water source, and do not expose the instrument directly to sunlight;
4. Do not place the instrument too close to a strong pyrogen or a strong radiation source (disinfection cabinet, induction cooker, microwave oven, etc.), which will affect the normal use of the instrument;
5. The instrument must use the original qualified accessories provided by the manufacturer;
6. After each operation, the operating handle should be cleaned and disinfected and kept in a safe place;
7. If you are not using the instrument, please turn off the power, dial out the instrument plug, and place all kinds of accessories.

## Packing list:

Main machine	1	Power cord	1
Big handle	1	Small handle	1
V10 handle	1		